

INTOXILYZER 8000 INSPECTION

Serial No.: 80-004945

Location: McIntosh Co

- Ashley
Check When Done
- A. Is the warm-up time less than 20 minutes? ✓
 - B. Is a three-pronged grounded outlet used? ✓
 - C. Is the breath tube heated? ✓
 - D. Is the diagnostic check complete? ✓
 - E. Is the time, date, and year correct (re-set if necessary)? ✓
 - F. Print test. (Attach test record.) ✓
 - G. Low AC. Use <0.03 AC in ACA mode. (Attach test record.)
Sim SN DR5132 Lot No. — Rep AC 0.021 AC ✓
 - H. Linearity Test. Use ≥ 0.25 AC in ACA mode. (Attach test record.)
Sim SN DR3375 Lot No. — Rep AC 0.292 AC ✓
 - I. Interferent check. Use 0.05% acetone plus 0.10 AC ethanol in ABA mode.
(Attach test record.)
Sim SN DR3454 Lot No. — Rep AC 0.10 AC + Acetone ✓
 - J. RFI check. Run CMS mode. Key radio on first room air.
(Attach test record.) ✓
 - K. Calibration Check. Use Ethanol Breath Standard cylinder. Do three sets in ACA mode.
(Attach test record.)

Lot No. 709352 Cyl No. 32 Exp. Date 12/23/2011

Test 1 <u>0.079</u>	Test 4 <u>0.079</u>	Test 7 <u>0.079</u>
Test 2 <u>0.079</u>	Test 5 <u>0.080</u>	Test 8 <u>0.079</u>
Test 3 <u>0.080</u>	Test 6 <u>0.080</u>	Test 9 <u>0.079 0.079</u>

Average 0.079

L. Tank pressure, Level 3 Func D Sub G P: 1057 psi Reg 1000 psi

M. Remarks/Maintenance record: Annual Inspection; Forms Load

2 June 2010
Date

Mary Plason
Field Inspector's Signature

14 June 2010
Date

Deb Kashner
Reviewed by

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004945
Location = ASHL 8164.13.00 06/09
06/02/2010 12:56

***** Printer Test *****

abcdefghijklmnopqrstuvwxyz1234567890-=_|
ABCDEFGHIJKLMNopqrstuvwxyz!@#%&^&*()_+?

abcdefghijklmnopqrstuvwxyz1234567890-=_|
ABCDEFGHIJKLMNopqrstuvwxyz!@#%&^&*()_+?

Current Instrument Setup


Data Entry Mode: Enabled
Start Test Sequence: DABACABA
Display Prelim Rslt? Yes
Display Third Digit? Yes
Inhib Printer(Y/N)? No
Display Volume? No
Disable On Memfull? No
of Print Copies? 1
Select Std (D/W/I)? Dry
Standard Value? 0.080
Standard Lot #? 709352
Standard Cyl #? 32
Standard Expiration? 12/23/2011
Oper No? 107501

Flow Cal. Date: 01/18/2010
Slope 709
Intercept -721404

IR Calibration Date: 01/18/2010
3um 9um

0th Coef(*100): -31758 -31251
1st Coef(*100): 257820 138163
2nd Coef(*100): 2248 1218
H2O adj(mg/l*10k): 753 539

***** Printer Test End *****

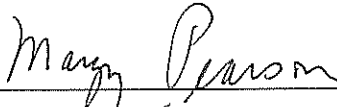

Operator Signature
MARGY PEARSON

Remarks:

Print

Form 106-I8000

80-004945
P. 39 of 46
15/2

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501


CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004945
Location = ASHL 8164.13.00 06/09
06/02/2010 13:19

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	13:20
02 Std. Sol.	0.016	13:23
03 Room Air	0.000	13:23
04 Std. Sol.	0.016	13:24
05 Room Air	0.000	13:24
06 Std. Sol.	0.017	13:25
07 Room Air	0.000	13:26

08 Sim Temp = 34.0°C

Simul Ser No = 000000
Std Sol No = 0000000
County = 26 Oper No. = 107501

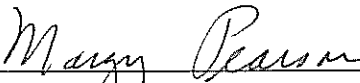

Operator Signature
MARGY PEARSON

Remarks:

Low AC 0.021 AC

Form 106-I8000

*80-004945
p. 46 of 46
JK*

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004945
Location = ASHL 8164.13.00 06/09
06/02/2010 13:12

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	13:13
02 Std. Sol.	0.292	13:14
03 Room Air	0.000	13:15
04 Std. Sol.	0.293	13:16
05 Room Air	0.000	13:16
06 Std. Sol.	0.293	13:17
07 Room Air	0.000	13:17


08 Sim Temp = 34.0°C

Simul Ser No = 000000

Std Sol No = 0000000

County = 26

Oper No. = 107501

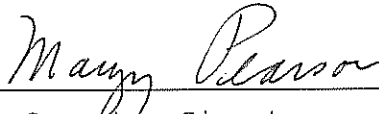

Operator Signature
MARGY PEARSON

Remarks:

Linearity 0.292 AC

Form 106-I8000

80-004945
P. 41 of 46
AM


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004945
Location = ASHL 8164.13.00 06/09
06/02/2010 13:26

Test	AC	Time
01 Room Air	0.000	13:27
02 *Subject Test	INT*	13:28
03 Room Air	0.000	13:28

*Invalid Test
Interferent Detected

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 11/01/1982
Sub Sex = Female Weight = 150
Test = DUI Cit = NA
Dr. Lic. = ND/DIS821456
Simul Ser No = 000000
Std Sol No = 0000000
County = 26 Oper No. = 107501



Operator Signature
MARGY PEARSON

Remarks:

Interferent Alc + Acetone

Form 106-I8000

80-004945
p. 42 of 46
JPM

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501


CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004945
Location = ASHL 8164.13.00 06/09
06/02/2010 13:33

Test	AC	Time
01 Diagnostic	OK	13:34
02 Room Air	RFI*	13:34
03 Room Air	0.000	13:35

*Invalid Test
 Inhibited - RFI

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 11/01/1982
Sub Sex = Female Weight = 150
Test = DUI Cit = NA
Dr. Lic. = ND/DIS821456
Lot No = 709352
Cyl No = 32
Expiration Date = 12/23/2011
County = 26 Oper No. = 107501

I followed the Approved Method and the instructions displayed by the Intoxilyzer in conducting this test.

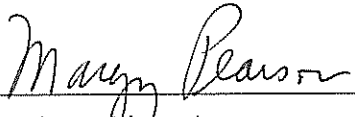

Operator Signature
MARGY PEARSON

Remarks:

RFI

Form 106-I8000

80-004945
p. 43 of 46
FM


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004945
Location = ASHL 8164.13.00 06/09
06/02/2010 13:30

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	13:30
02 Std. Gas	0.079	13:30
03 Room Air	0.000	13:31
04 Std. Gas	0.079	13:31
05 Room Air	0.000	13:32
06 Std. Gas	0.080	13:32
07 Room Air	0.000	13:33

Lot No = 709352
Cyl No = 32
Exp Date = 12/23/2011
County = 26 Oper No. = 107501



Operator/Signature
MARGY PEARSON

Remarks:

Cal Check AEA 0.080 AC

Form 106-I8000

*80-004945
P. 444 46
FM*


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004945
Location = ASHL 8164.13.00 06/09
06/02/2010 13:35

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	13:36
02 Std. Gas	0.079	13:36
03 Room Air	0.000	13:37
04 Std. Gas	0.080	13:37
05 Room Air	0.000	13:38
06 Std. Gas	0.080	13:38
07 Room Air	0.000	13:38

Lot No = 709352
Cyl No = 32
Exp Date = 12/23/2011
County = 26 Oper No. = 107501

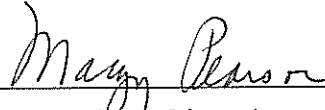

Operator Signature
MARGY PEARSON

Remarks:

Cal Check ACA

Form 106-I8000

80-004945
P. 45 of 46
KM


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004945
Location = ASHL 8164.13.00 06/09
06/02/2010 13:39

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	13:40
02 Std. Gas	0.079	13:40
03 Room Air	0.000	13:41
04 Std. Gas	0.079	13:41
05 Room Air	0.000	13:42
06 Std. Gas	0.079	13:42
07 Room Air	0.000	13:43

Lot No = 709352
Cyl No = 32
Exp Date = 12/23/2011
County = 26 Oper No. = 107501

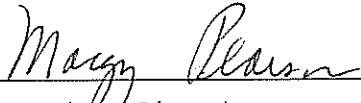

Operator Signature
MARGY PEARSON

Remarks:

Cal Check AC A 0.080 AC

Form 106-I8000

80-004945
P. 96 of 46
JEM