



INTOXILYZER® 8000 INSTALLATION AND REPAIR CHECKOUT
 OFFICE OF ATTORNEY GENERAL
 CRIME LABORATORY DIVISION
 Toxicology Section/Breath Alcohol Program
 SFN 59281 (1-2010)

Serial Number: <u>80-004952</u>	Instrument Location: ^{MRP} <u>± Ransom Co SO LISBON</u>
Reason for Install/Repair:	
<input checked="" type="checkbox"/> Install After Lab Repair/ <u>Inspection</u> <input type="checkbox"/> Other (Specify) _____	

Check When Done:

- A. Surge Protector Installed/Properly Grounded.
- B. Telephone Line Connected to Intoxilyzer® 8000.
- C. Breath Tube Heated.
- D. Review/Enter Preliminary Data Entry (Level 2, Function E).
- E. Review/Enter Gas Setup (Level 1, Function S).
- F. Print Test (Level 1, Function P). (Attach Test Record.)
- G. ACA Test (Level 1, Function C). (Attach Test Record.)
- H. RFI Test (Use CMS Mode-First Room Air). (Attach Test Record.)
- I. Repair and/or Maintenance Performed: Install - Swapped out Intox 5000 G

Margy Pearson
 Field Inspector Signature

2 June 2010
 Date

Deb Kashner
 Reviewed By

14 June 2010
 Date

80-004952
 P. 36 of 39
 km

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004952
Location = LISB 8164.13.00 06/09
06/02/2010 17:40

***** Printer Test *****

abcdefghijklmnopqrstuvwxy1234567890-|=|
ABCDEFGHIJKLMN0PQRSTUVWXYZ!@#\$\$%^&*()_+?

abcdefghijklmnopqrstuvwxy1234567890-|=|
ABCDEFGHIJKLMN0PQRSTUVWXYZ!@#\$\$%^&*()_+?

Current Instrument Setup


Data Entry Mode: Enabled
Start Test Sequence: DABACABA
Display Prelim Rslt? Yes
Display Third Digit? Yes
Inhib Printer(Y/N)? No
Display Volume? No
Disable On Memfull? Yes
of Print Copies? 1
Select Std (D/W/I)? Dry
Standard Value? 0.080
Standard Lot #? 709352
Standard Cyl #? 22
Standard Expiration? 12/23/2011
Oper No? 107501

Flow Cal. Date: 01/19/2010
Slope 694
Intercept -638616

IR Calibration Date: 01/19/2010
 3um 9um

0th Coef(*100): -32832 -34603
1st Coef(*100): 262566 131011
2nd Coef(*100): 1418 975
H2O adj(mg/l*10k): 686 550

***** Printer Test End *****

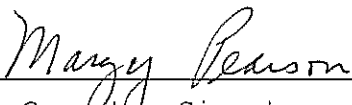

Operator Signature
MARGY PEARSON

Remarks:

Print


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004952
Location = LISB 8164.13.00 06/09
06/02/2010 17:52

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	17:53
02 Std. Gas	0.079	17:53
03 Room Air	0.000	17:53
04 Std. Gas	0.079	17:54
05 Room Air	0.000	17:54
06 Std. Gas	0.079	17:55
07 Room Air	0.000	17:55

Lot No = 709352
Cyl No = 22
Exp Date = 12/23/2011
County = 37 Oper No. = 107501

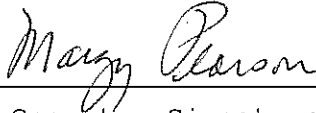

Operator Signature
MARGY PEARSON

Remarks:

ACA I8000 Install

Form 106-I8000

80-004952
P. 38 of 39 KM

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004952
Location = LISB 8164.13.00 06/09
06/02/2010 17:56

Test	AC	Time
01 Diagnostic	OK	17:57
02 Room Air	RFI*	17:57
03 Room Air	0.000	17:58

*Invalid Test
Inhibited - RFI

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 11/01/1982
Sub Sex = Female Weight = 150
Test = DUI Cit = NA
Dr. Lic. = ND/DIS821456
Lot No = 709352
Cyl No = 22
Expiration Date = 12/23/2011
County = 37 Oper No. = 107501

I followed the Approved Method and the instructions displayed by the Intoxilyzer in conducting this test.



Operator Signature
MARGY PEARSON

Remarks:

RFI Test

Form 106-I8000

80-004952
P.39 of 39
JAM