

INTOXILYZER 8000 INSPECTION

Serial No.: 80-004956

Location: Emmons Co SO

Lint on
Check When Done

- A. Is the warm-up time less than 20 minutes? ✓
- B. Is a three-pronged grounded outlet used? ✓
- C. Is the breath tube heated? ✓
- D. Is the diagnostic check complete? ✓
- E. Is the time, date, and year correct (re-set if necessary)? ✓
- F. Print test. (Attach test record.) ✓
- G. Low AC. Use <0.03 AC in ACA mode. (Attach test record.)
Sim SN DR5132 Lot No. — Rep AC 0.021 AC ✓
- H. Linearity Test. Use ≥ 0.25 AC in ACA mode. (Attach test record.)
Sim SN DR3375 Lot No. — Rep AC 0.292 AC ✓
- I. Interferent check. Use 0.05% acetone plus 0.10 AC ethanol in ABA mode.
(Attach test record.)
Sim SN DR3454 Lot No. — Rep AC 0.10 AC + Acetone ✓
- J. RFI check. Run CMS mode. Key radio on first room air.
(Attach test record.) ✓
- K. Calibration Check. Use Ethanol Breath Standard cylinder. Do three sets in ACA mode.
(Attach test record.)

Lot No. 709352 Cyl No. 16 Exp. Date 12/23/11

Test 1 0.077 Test 4 0.078 Test 7 0.077

Test 2 0.077 Test 5 0.077 Test 8 0.077

Test 3 0.077 Test 6 0.077 Test 9 0.077

Average 0.077 AC

L. Tank pressure, Level 3 Func D Sub G P: 770 psi Reg 780 psi

M. Remarks/Maintenance record: Annual Inspection MODEM/
Forms Load FAX 254 5311

2 June 2011
Date

Mary Pearson
Field Inspector's Signature

14 June 2010
Date

Deb Keshner
Reviewed by

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004956
Location = LINT 8164.13.00 06/09
06/02/2010 10:28

***** Printer Test *****

abcdefghijklmnopqrstuvwxy1234567890-==|
ABCDEFGHIJKLMNPOQRSTUVWXYZ!@#%&^*()_+?

abcdefghijklmnopqrstuvwxy1234567890-==|
ABCDEFGHIJKLMNPOQRSTUVWXYZ!@#%&^*()_+?

Current Instrument Setup


Data Entry Mode: Enabled
Start Test Sequence: DABACABA
Display Prelim Rslt? Yes
Display Third Digit? Yes
Inhib Printer(Y/N)? No
Display Volume? No
Disable On Memfull? No
of Print Copies? 1
Select Std (D/W/I)? Dry
Standard Value? 0.080
Standard Lot #? 709352
Standard Cyl #? 16
Standard Expiration? 12/23/2011
Oper No? 107501

Flow Cal. Date: 01/25/2010
Slope 672
Intercept -627011

IR Calibration Date: 02/12/2010
 3um 9um

0th Coef(*100): -309 61
1st Coef(*100): 241249 120296
2nd Coef(*100): 3156 1700
H2O adj(mg/l*10k): 555 539

***** Printer Test End *****



Operator Signature
MARGY PEARSON

Remarks:

Print Test

Form 106-I8000

80-004956
P. 41 of 48
JW

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004956
Location = LINT 8164.13.00 06/09
06/02/2010 10:39

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	10:40
02 Std. Sol.	0.020	10:41
03 Room Air	0.000	10:41
04 Std. Sol.	0.020	10:42
05 Room Air	0.000	10:42
06 Std. Sol.	0.021	10:43
07 Room Air	0.000	10:44


08 Sim Temp = 34.0°C

Simul Ser No = 0000000

Std Sol No = 000000

County = 15

Oper No. = 107501

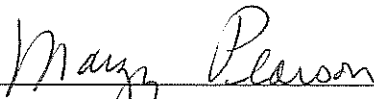

Operator Signature
MARGY PEARSON

Remarks:

Low AC 0.021 AC

Form 106-I8000

80-004956
P. 42 of 48
PM

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004956
Location = LINT 8164.13.00 06/09
06/02/2010 10:33

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	10:34
02 Std. Sol.	0.283	10:35
03 Room Air	0.000	10:35
04 Std. Sol.	0.285	10:36
05 Room Air	0.000	10:37
06 Std. Sol.	0.285	10:37
07 Room Air	0.000	10:38


08 Sim Temp = 34.0°C

Simul Ser No = 0000000

Std Sol No = 000000

County = 15

Oper No. = 107501

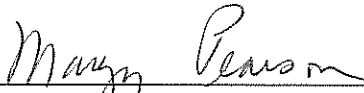

Operator Signature
MARGY PEARSON

Remarks:

Linearity 0.292 AC

Form 106-I8000

80-004956
p. 43 of 48
AM


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004956
Location = LINT 8164.13.00 06/09
06/02/2010 10:47

Test	AC	Time
01 Room Air	0.000	10:48
02 *Subject Test	INT*	10:48
03 Room Air	0.000	10:49

*Invalid Test
Interferent Detected

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 11/01/1982
Sub Sex = Female Weight = 150
Test = DUI Cit = NA
Dr. Lic. = ND/DIS821456
Simul Ser No = 0000000
Std Sol No = 000000
County = 15 Oper No. = 107501

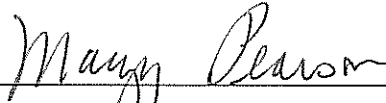

Operator Signature
MARGY PEARSON

Remarks:

Interferent - Alc + Acetone (0.05%)

Form 106-I8000

80-004956
P. 44 of 48
PM

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501


CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004956
Location = LINT 8164.13.00 06/09
06/02/2010 10:49

Test	AC	Time
01 Diagnostic	OK	10:50
02 Room Air	RFI*	10:51
03 Room Air	0.000	10:51

*Invalid Test
Inhibited - RFI

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 11/01/1982
Sub Sex = Female Weight = 150
Test = DUI Cit = NA
Dr. Lic. = ND/DIS821456
Simul Ser No = 0000000
Std Sol No = 000000
County = 15 Oper No. = 107501

I followed the Approved Method and the instructions displayed by the Intoxilyzer in conducting this test.



Operator Signature
MARGY PEARSON

Remarks:

RFI

Form 106-I8000

80-004956
P. 45 of 48
JAM


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004956
Location = LINT 8164.13.00 06/09
06/02/2010 10:52

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	10:53
02 Std. Gas	0.077	10:53
03 Room Air	0.000	10:54
04 Std. Gas	0.077	10:54
05 Room Air	0.000	10:55
06 Std. Gas	0.077	10:55
07 Room Air	0.000	10:56

Lot No = 709352
Cyl No = 16
Exp Date = 12/23/2011
County = 15 Oper No. = 107501



Operator Signature
MARGY PEARSON

Remarks:

Cal Check 0.080 AC

Form 106-I8000

*80-004956
P. 46448
FM*

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004956
Location = LINT 8164.13.00 06/09
06/02/2010 10:56

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	10:57
02 Std. Gas	0.078	10:57
03 Room Air	0.000	10:58
04 Std. Gas	0.077	10:58
05 Room Air	0.000	10:59
06 Std. Gas	0.077	10:59
07 Room Air	0.000	10:59

Lot No = 709352
Cyl No = 16
Exp Date = 12/23/2011
County = 15 Oper No. = 107501

Margy Pearson

Operator Signature
MARGY PEARSON

Remarks:

Cal Check 0.080 AC

Form 106-18000

80-004956
P. 47 of 48
JEM


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004956
Location = LINT 8164.13.00 06/09
06/02/2010 11:00

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	11:01
02 Std. Gas	0.077	11:01
03 Room Air	0.000	11:01
04 Std. Gas	0.077	11:02
05 Room Air	0.000	11:02
06 Std. Gas	0.077	11:03
07 Room Air	0.000	11:03

Lot No = 709352
Cyl No = 16
Exp Date = 12/23/2011
County = 15 Oper No. = 107501

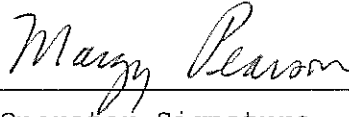

Operator Signature
MARGY PEARSON

Remarks:

Cal Check 0.080AC

Form 106-I8000

80-004956
P. 48 of 48
JAN