

INTOXILYZER® 8000 INSPECTION

Serial No.: 80-002667

Location: LANG (CAVALIER COUNTY S.O.)

- Check When Done
- A. Is the breath tube heated? ✓
 - B. Is the diagnostic check complete? ✓
 - C. Is the time, date, and year correct (re-set if necessary)? ✓
 - D. Print test. (Attach test record.) ✓
 - E. Low AC. Use <0.03 AC in ACA mode. (Attach test record.)
Sim SN DR5115 Lot No. 1211 Rep AC 0.020 ✓
 - F. High AC. Use ≥ 0.25 AC in ACA mode. (Attach test record.)
Sim SN DR5143 Lot No. 13080 Rep AC 0.300 ✓
 - G. Interferent check. Use 0.05% acetone plus 0.10 AC ethanol in ABA mode.
(Attach test record.)
Sim SN DR3375 Lot No. IGS#2 Rep AC 0.100AC + 0.05% ACETONE ✓
 - H. RFI check. Run CMS mode. Key radio on first room air.
(Attach test record.) ✓

I. Calibration Check. Use Ethanol Breath Standard cylinder. Do three sets in ACA mode.
(Attach test record.)

Lot No. <u>16512080A2</u>	Cyl No. <u>060</u>	Exp. Date <u>08/01/2014</u>
Test 1 <u>0.080</u>	Test 4 <u>0.080</u>	Test 7 <u>0.079</u>
Test 2 <u>0.079</u>	Test 5 <u>0.080</u>	Test 8 <u>0.080</u>
Test 3 <u>0.080</u>	Test 6 <u>0.079</u>	Test 9 <u>0.080</u>
Average <u>0.080</u>		

J. Tank pressure, Level 3 Func D Sub G P: 391 psi Reg 400 psi

K. Remarks/Maintenance record: REPLACED BREATH HOSE O-RING, REPLACED SIMULATOR RETURN O-RING, ANNUAL INSPECTION

02 JUN 14

Date

Charles E. Ed

Field Inspector's Signature

04 June 14 ✓

Date

Deb Kashner

Reviewed by

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-002667
Location = LANG 8164.13.00 06/09
06/02/2014 11:18

***** Printer Test *****

abcdefghijklmnopqrstuvwxyz1234567890-=_|
ABCDEFGHIJKLMNopqrstuvwxyz!@#\$\$%^&*()_+?

abcdefghijklmnopqrstuvwxyz1234567890-=_|
ABCDEFGHIJKLMNopqrstuvwxyz!@#\$\$%^&*()_+?

Current Instrument Setup

Data Entry Mode: Enabled
Start Test Sequence: DABACABA
Display Prelim Rslt? Yes
Display Third Digit? Yes
Inhib Printer(Y/N)? No
Display Volume? No
Disable On Memfull? Yes
of Print Copies? 1
Select Std (D/W/I)? Dry
Standard Value? 0.080
Standard Lot #? 16512080A2
Standard Cyl #? 60
Standard Expiration? 08/01/2014
Oper No? 666666

Flow Cal. Date: 09/10/2007
Slope 693
Intercept -570251

IR Calibration Date: 06/03/2011
 3um 9um

0th Coef(*100): -3044 -42
1st Coef(*100): 248787 127659
2nd Coef(*100): 6023 1886
H2O adj(mg/l*10k): 632 468

***** Printer Test End *****



Operator Signature
CHARLES EDER

Remarks:

PRINT TEST

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-002667
Location = LANG 8164.13.00 06/09
06/02/2014 11:36

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	11:37
02 Std. Sol.	0.017	11:37
03 Room Air	0.000	11:38
04 Std. Sol.	0.019	11:39
05 Room Air	0.000	11:39
06 Std. Sol.	0.019	11:40
07 Room Air	0.000	11:40


08 Sim Temp = 34.0°C

Simul Ser No = DR5115

Std Sol No = 12111

County = 10

Oper No. = 666666

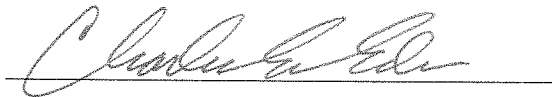

Operator Signature

CHARLES EDER

Remarks:

LOW AC
0.020 AC

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501


CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-002667
Location = LANG 8164.13.00 06/09
06/02/2014 11:42

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	11:42
02 Std. Sol.	0.291	11:43
03 Room Air	0.000	11:44
04 Std. Sol.	0.300	11:44
05 Room Air	0.000	11:45
06 Std. Sol.	0.305	11:46
07 Room Air	0.000	11:46

08 Sim Temp = 34.0°C

Simul Ser No = DR5143
Std Sol No = 13080
County = 10 Oper No. = 666666

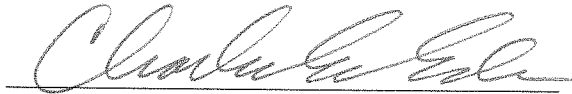

Operator Signature
CHARLES EDER

Remarks:

HIGH AC
0.300 AC

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-002667
Location = LANG 8164.13.00 06/09
06/02/2014 11:32

Test	AC	Time
01 Room Air	0.000	11:33
02 *Subject Test	INT*	11:34
03 Room Air	0.000	11:35

*Invalid Test
Interferent Detected

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 01/01/1982
Sub Sex = Female Weight = 150
Test = DUI Cit = INTERFERENT
Dr. Lic. = ND/DIS821456
Lot No = 16512080A2
Cyl No = 60
Expiration Date = 08/01/2014
County = 10 Oper No. = 666666



Operator Signature
CHARLES EDER

Remarks:

*INTERFERENT
CHECK*

0.100AC + 0.05% ACETONE

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501


CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-002667
Location = LANG 8164.13.00 06/09
06/02/2014 11:21

Test	AC	Time
01 Diagnostic	OK	11:23
02 Room Air	RFI*	11:23
03 Room Air	0.000	11:23

*Invalid Test
Inhibited - RFI

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 01/01/1982
Sub Sex = Female Weight = 150
Test = DUI Cit = RFI CHECK
Dr. Lic. = ND/DIS821456
Lot No = 16512080A2
Cyl No = 60
Expiration Date = 08/01/2014
County = 10 Oper No. = 666666

I followed the Approved Method and the instructions displayed by the Intoxilyzer in conducting this test.

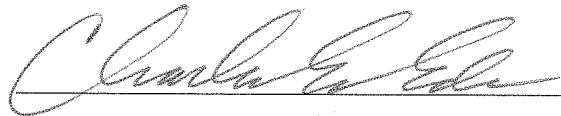

Operator Signature
CHARLES EDER

Remarks:

RFI CHECK

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-002667
Location = LANG 8164.13.00 06/09
06/02/2014 11:24

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	11:25
02 Std. Gas	0.080	11:25
03 Room Air	0.000	11:26
04 Std. Gas	0.079	11:26
05 Room Air	0.000	11:26
06 Std. Gas	0.080	11:27
07 Room Air	0.000	11:27

Lot No = 16512080A2
Cyl No = 60
Exp Date = 08/01/2014
County = 10 Oper No. = 666666



Operator Signature
CHARLES EDER

Remarks: *CALIBRATION*
Check
0.080AC
Form 106-I8000


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-002667
Location = LANG 8164.13.00 06/09
06/02/2014 11:28

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	11:29
02 Std. Gas	0.080	11:29
03 Room Air	0.000	11:30
04 Std. Gas	0.080	11:30
05 Room Air	0.000	11:31
06 Std. Gas	0.079	11:31
07 Room Air	0.000	11:31

Lot No = 16512080A2
Cyl No = 60
Exp Date = 08/01/2014
County = 10 Oper No. = 666666

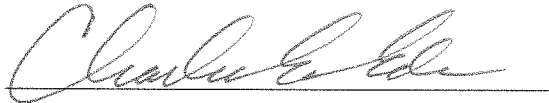

Operator Signature
CHARLES EDER

Remarks: *CALIBRATION
CHECK
0.080 AC*
Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-002667
Location = LANG 8164.13.00 06/09
06/02/2014 11:47

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	11:48
02 Std. Gas	0.079	11:48
03 Room Air	0.000	11:49
04 Std. Gas	0.080	11:49
05 Room Air	0.000	11:50
06 Std. Gas	0.080	11:50
07 Room Air	0.000	11:50

Lot No = 16512080A2
Cyl No = 60
Exp Date = 08/01/2014
County = 10 Oper No. = 666666



Operator Signature
CHARLES EDER

Remarks:

CALIBRATION
CHECK
0.080 AC

Form 106-I8000