



**ETHANOL GAS STANDARD CYLINDER REPORT**  
OFFICE OF ATTORNEY GENERAL  
CRIME LABORATORY DIVISION  
Toxicology Section/Breath Alcohol Program  
SFN 59282 (10-2013)

Chemical Test Operator Name (Print) <b>DEUBERT HOFFMANN</b>		
Location <b>SHER</b>	Intoxilyzer® Serial Number <b>80-005950</b>	
Gas Lot Number <b>09216080A3</b>	Gas Cylinder Number <b>30</b>	Gas Expiration Date <b>05-05-2018</b>

Check When Done:

- 1. Scan/Enter Gas Cylinder Information (Level 1, Function S).
- 2. Perform an ACA Test (Level 1, Function C).  
Write in the ACA Test Results (Reported to 3 Digits, Ex. 0.081).
  - A. **0.082** AC
  - B. **0.082** AC
  - C. **0.082** AC
- 3. Complete the Top Portion of the Intoxilyzer® Record (SFN50496, Form 120-G) and Place it by the Intoxilyzer® for Use.
- 4. File Previous Intoxilyzer® Record (SFN50496, Form 120-G) at the Intoxilyzer® Location at the Agency.
- 5. Send the Following to the Crime Laboratory:
  - A. Ethanol Gas Standard Cylinder Report (SFN59282, Form 105-G).
  - B. ACA Test Record.

Chemical Test Operator Signature <b>Debert Hoffmann</b>	Date <b>06-29-2016 ✓</b>
Reviewed By (Crime Laboratory Use Only) <b>Deb Shanaver</b>	Date <b>06 July 16</b>

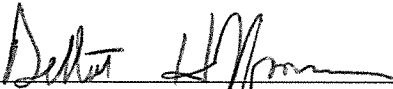
Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-005950  
Location = SHER      8164.13.00 06/09  
06/29/2016      17:20

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	17:21
02 Std. Gas	0.082	17:21
03 Room Air	0.000	17:22
04 Std. Gas	0.082	17:22
05 Room Air	0.000	17:22
06 Std. Gas	0.082	17:23
07 Room Air	0.000	17:23

Lot No = 09216080A3  
Cyl No = 30  
Exp Date = 05/05/2018  
County = 42      Oper No. = 048909

  
\_\_\_\_\_  
Operator Signature  
DELBERT HOFFMANN

Remarks:      *NCA*