

INTOXILYZER 8000 INSPECTION

Serial No.: 80-004188

Location: TOXL

- | | Check When Done |
|--|-----------------|
| A. Is the warm-up time less than 20 minutes? | <u>✓</u> |
| B. Is a three-pronged grounded outlet used? | <u>✓</u> |
| C. Is the breath tube heated? | <u>✓</u> |
| D. Is the diagnostic check complete? | <u>✓</u> |
| E. Is the time, date, and year correct (re-set if necessary)? | <u>✓</u> |
| F. Print test. (Attach test record.) | <u>✓</u> |
| G. Low AC. Use <0.03 AC in ACA mode. (Attach test record.)
Sim SN <u>DR5132</u> Lot No. <u> </u> Rep AC <u>0.021</u> | <u>✓</u> |
| H. Linearity Test. Use ≥ 0.25 AC in ACA mode. (Attach test record.)
Sim SN <u>DR5185</u> Lot No. <u> </u> Rep AC <u>0.297</u> | <u>✓</u> |
| I. Interferent check. Use 0.05% acetone plus 0.10 AC ethanol in ABA mode.
(Attach test record.)
Sim SN <u>DR3846</u> Lot No. <u> </u> Rep AC <u>0.10 AC + acetone</u> | <u>✓</u> |
| J. RFI check. Run CMS mode. Key radio on first room air.
(Attach test record.) | <u>✓</u> |
| K. Calibration Check. Use Ethanol Breath Standard cylinder. Do three sets in ACA mode.
(Attach test record.) | |

Lot No. <u>1541008041</u>	Cyl No. <u>25</u>	Exp. Date <u>07/01/2012</u>
Test 1 <u>0.079</u>	Test 4 <u>0.079</u>	Test 7 <u>0.080</u>
Test 2 <u>0.079</u>	Test 5 <u>0.079</u>	Test 8 <u>0.080</u>
Test 3 <u>0.079</u>	Test 6 <u>0.080</u>	Test 9 <u>0.079</u>
Average <u>0.079</u>		

L. Tank pressure, Level 3 Func D Sub G P: 391 psi Reg 400 psi

M. Remarks/Maintenance record: Annual inspection after factory EEPROM change.

25 July 2011
Date
7-26-11
Date

Deb Kashner
Field Inspector's Signature
[Signature]
Reviewed by

80-004188
p. 152 of 160
SM

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004188
Location = TOXL 8164.13.00 06/09
07/25/2011 09:59

***** Printer Test *****

abcdefghijklmnopqrstuvwxy1234567890-=_|
ABCDEFGHIJKLMNPOQRSTUVWXYZ!@#\$\$%^&*()_+?

abcdefghijklmnopqrstuvwxy1234567890-=_|
ABCDEFGHIJKLMNPOQRSTUVWXYZ!@#\$\$%^&*()_+?

Current Instrument Setup

Data Entry Mode: Enabled
Start Test Sequence: DABACABA
Display Prelim Rslt? Yes
Display Third Digit? Yes
Inhib Printer(Y/N)? No
Display Volume? No
Disable On Memfull? Yes
of Print Copies? 1
Select Std (D/W/I)? Dry
Standard Value? 0.080
Standard Lot #? 15410080A1
Standard Cyl #? 25
Standard Expiration? 07/01/2012
Oper No? 777777

Flow Cal. Date: 07/05/2011
Slope 668
Intercept -597204

IR Calibration Date: 07/05/2011
 3um 9um

0th Coef(*100): -19325 -21653
1st Coef(*100): 277687 137181
2nd Coef(*100): 2891 1139
H2O adj(mg/l*10k): 555 619

***** Printer Test End *****

Deb Kasher

Operator Signature
DEB KASHUR

Remarks:

Print test

Form 106-I8000

80-004188
P.153 of 160
KTM

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004188
Location = TOXL 8164.13.00 06/09
07/25/2011 10:01

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	10:02
02 Std. Sol.	0.020	10:02
03 Room Air	0.000	10:03
04 Std. Sol.	0.020	10:04
05 Room Air	0.000	10:04
06 Std. Sol.	0.021	10:05
07 Room Air	0.000	10:05

08 Sim Temp = 34.0°C

Simul Ser No = 0

Std Sol No = 0

County = 08

Oper No. = 777777

Deb Kashur
Operator Signature
DEB KASHUR

Remarks:

low AC 0.021 AC

Form 106-I8000

80-004188
P.154 of 160
PM

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004188
Location = TOXL 8164.13.00 06/09
07/25/2011 10:56

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	10:56
02 Std. Sol.	0.287	10:57
03 Room Air	0.000	10:58
04 Std. Sol.	0.288	10:58
05 Room Air	0.000	10:59
06 Std. Sol.	0.289	11:00
07 Room Air	0.000	11:00

08 Sim Temp = 34.0°C

Simul Ser No = 0

Std Sol No = 0

County = 08

Oper No. = 777777

Deb Kasher
Operator Signature
DEB KASHUR

Remarks:

Linearity test 0.297 AC

Form 106-I8000

*80-004188
p. 155 of 160
FM*

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004188
Location = TOXL 8164.13.00 06/09
07/25/2011 17:01

Test	AC	Time
01 Room Air	0.000	17:01
02 *Subject Test	INT*	17:02
03 Room Air	0.000	17:02

*Invalid Test
Interferent Detected

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 02/01/1992
Sub Sex = Female Weight = 150
Test = DUI Cit = NA
Dr. Lic. = ND/DIS921456
Lot No = 15410080A1
Cyl No = 25
Expiration Date = 07/01/2012
County = 08 Oper No. = 777777

Deb Kashur

Operator Signature
DEB KASHUR

Remarks:

Interferent check 0.10Ac + acetone

Form 106-I8000

80-004188
P. 156 of 160
AM

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004188
Location = TOXL 8164.13.00 06/09
07/25/2011 17:03

Test	AC	Time
01 Diagnostic	OK	17:04
02 Room Air	RFI*	17:04
03 Room Air	0.000	17:04

*Invalid Test
Inhibited - RFI

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 02/01/1992
Sub Sex = Female Weight = 150
Test = DUI Cit = NA
Dr. Lic. = ND/DIS921456
Lot No = 15410080A1
Cyl No = 25
Expiration Date = 07/01/2012
County = 08 Oper No. = 777777

I followed the Approved Method and the instructions displayed by the Intoxilyzer in conducting this test.

Deb Kasher

Operator Signature
DEB KASHUR

Remarks:

RFI test

Form 106-I8000

80-004188
P 157A/60
PTM

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004188
Location = TOXL 8164.13.00 06/09
07/22/2011 11:22

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	11:22
02 Std. Gas	0.079	11:23
03 Room Air	0.000	11:23
04 Std. Gas	0.079	11:24
05 Room Air	0.000	11:24
06 Std. Gas	0.079	11:25
07 Room Air	0.000	11:25

Lot No = 15410080A1
Cyl No = 25
Exp Date = 07/01/2012
County = 08 Oper No. = 777777

Deb Kasher

Operator Signature
DEB KASHUR

Remarks:

Calibration check 0.080 AC

Form 106-I8000

*80-004188
0.158-160
DM*


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004188
Location = TOXL 8164.13.00 06/09
07/25/2011 09:51

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	09:52
02 Std. Gas	0.079	09:52
03 Room Air	0.000	09:52
04 Std. Gas	0.079	09:53
05 Room Air	0.000	09:53
06 Std. Gas	0.080	09:54
07 Room Air	0.000	09:54

Lot No = 15410080A1
Cyl No = 25
Exp Date = 07/01/2012
County = 08 Oper No. = 777777

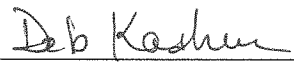

Operator Signature
DEB KASHUR

Remarks:
Calibration check 0.080 AC

Form 106-I8000

*80-004188
P.1594160
JEM*

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004188
Location = TOXL 8164.13.00 06/09
07/25/2011 09:55

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	09:56
02 Std. Gas	0.080	09:56
03 Room Air	0.000	09:56
04 Std. Gas	0.080	09:57
05 Room Air	0.000	09:57
06 Std. Gas	0.079	09:58
07 Room Air	0.000	09:58

Lot No = 15410080A1
Cyl No = 25
Exp Date = 07/01/2012
County = 08 Oper No. = 777777

Deb Kashur
Operator Signature
DEB KASHUR

Remarks:

Calibration check 0.080AC

Form 106-I8000

80-004188
P. 160 of 160
JEM