

OFFICE OF ATTORNEY GENERAL
 CRIME LABORATORY DIVISION
 Toxicology Section/Breath Alcohol Program
 Instrument Certification and Maintenance Unit

BrF-005-8

J Q. Flow Monitor. Level 3 Func D Sub F. R 170 F 0.000 L/s

✓ R. Temperature Monitor. Level 3 Func D Sub T
 C 47.0 A 49.0 I 0.0
 B 45.0 V 0.0

✓ S. Gas Tank Monitor. Level 3 Func D Sub G
✓ T. Configure Wet Bath Sim. Level 1 Func S
979 psi Intox 950 psi Reg

✓ U. Zero Alcohol Test Level 1 Func C
610329 Sim SN — Lot 0.000 Rep Value
0.000 AC 0.000 AC 0.000 AC

✓ V. Low AC Test. Level 1 Func C
DR519 Sim SN — Lot 0.021 Rep Value
0.017 AC 0.019 AC 0.020 AC

✓ W. Linearity Test. Level 1 Func C
DR513 Sim SN — Lot 0.297 Rep Value
0.289 AC 0.291 AC 0.292 AC

✓ X. Interferent Test.
 Level 1 Func B
 Acetone DR3845 Sim SN ✓ Lot 0.050% Conc
✓ Display ✓ Printout
 Methanol DR3846 Sim SN ✓ Lot 0.10% Conc
✓ Display ✓ Printout
 IPA DR5191 Sim SN ✓ Lot 0.10% Conc
✓ Display ✓ Printout

✓ Y. Ambient Test
✓ Display ✓ Printout

✓ Z. Invalid Sample Test
✓ Display ✓ Printout

✓ AA. RFI Detector Test. CMS Mode
✓ Display ✓ Printout

✓ BB. Configure Dry Gas Standard. Level 1 Func S

✓ CC. Stability Test at 0.080AC. Level 3 Func T Sub S
0.0809 Ave. 0.0005 Std Dev 0.619% Rel Std Dev

DD. Maintenance/Repairs

26 August 2009
 Date
04 Sept 2009
 Date

Mary Plam
 Signature
Deb Kashner
 Reviewed by

*P. 2 of 176
 FM*

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-002669
Location = TOXL 8164.12.00 06/09
08/26/2009 08:27

***** Printer Test *****

abcdefghijklmnopqrstuvwxyz1234567890-=
ABCDEFGHIJKLMNopqrstuvwxyz!@#%&*()_+?

abcdefghijklmnopqrstuvwxyz1234567890-=
ABCDEFGHIJKLMNopqrstuvwxyz!@#%&*()_+?

Current Instrument Setup


Data Entry Mode: Enabled
Start Test Sequence: DABACABA
Display Prelim Rslt? Yes
Display Third Digit? Yes
Inhib Printer(Y/N)? No
Display Volume? No
Disable On Memfull? Yes
of Print Copies? 1
Select Std (D/W/I)? Wet
Standard Value? 0.110
Standard Solution #? DR00000
Sim Serial #? 1
Oper No? 107501

Flow Cal. Date: 09/10/2007
Slope 702
Intercept -608746

IR Calibration Date: 09/10/2007
 3um 9um

0th Coef(*100): -4604 -8890
1st Coef(*100): 254947 133181
2nd Coef(*100): 4056 1165
H2O adj(mg/l*10k): 616 347

***** Printer Test End *****

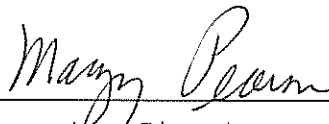

Operator Signature
MARGY PEARSON

Remarks:

Print Test


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-002669
Location = TOXL 8164.12.00 06/09
08/26/2009 08:59

Test	AC	Time
01 Room Air	0.000	09:01
02 *Subject Test	REF*	09:01
03 Room Air	0.000	09:01

*Subject Test Refused

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 11/01/1982
Sub Sex = Female Weight = 150
Test = DUI Cit = NA
Dr. Lic. = ND/DIS821456
Simul Ser No = 1
Std Sol No = DR00000
County = 08 Oper No. = 107501

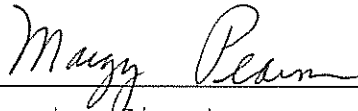

Operator Signature
MARGY PEARSON

Remarks:

Refusal

Form 106-I8000

*P. 4 of 16
JRM*


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-002669
Location = TOXL 8164.12.00 06/09
08/26/2009 09:02

DIAGNOSTIC

Voltage/Current Test	Pass
RAM Test	Pass
EEPROM Checksum Test	Pass
Real Time Clock Test	Pass
DSP Test	Pass
Analytical Stability Test	Pass
Modem Test	Pass
Temp Regulation Test	Pass

Oper No. = 107501

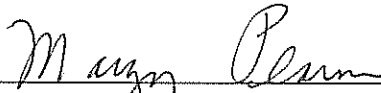

Operator Signature
MARGY PEARSON

Remarks:

DIAG

Form 106-I8000

P. Sep 16
JRM

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-002669
Location = TOXL 8164.12.00 06/09
08/26/2009 13:32

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	13:33
02 Std. Sol.	0.000	13:34
03 Room Air	0.000	13:34
04 Std. Sol.	0.000	13:35
05 Room Air	0.000	13:35
06 Std. Sol.	0.000	13:36
07 Room Air	0.000	13:37


08 Sim Temp = 34.0°C

Simul Ser No = 1

Std Sol No = DR00000

County = 08

Oper No. = 107501

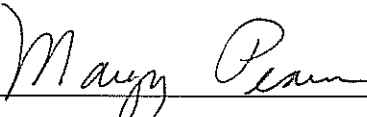

Operator Signature
MARGY PEARSON

Remarks:

Water 0.00Ac

Form 106-I8000

P. 2 of 16
JFM

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-002669
Location = TOXL 8164.12.00 06/09
08/26/2009 11:45

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	11:46
02 Std. Sol.	0.017	11:50
03 Room Air	0.000	11:50
04 Std. Sol.	0.019	11:51
05 Room Air	0.000	11:51
06 Std. Sol.	0.020	11:52
07 Room Air	0.000	11:53

08 Sim Temp = 34.0°C

Simul Ser No = 1
Std Sol No = DR00000
County = 08

Oper No. = 107501



Operator Signature
MARGY PEARSON

Remarks:

LOW 0.021 AC

Form 106-I8000

P. 7 of 16
[Signature]

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-002669
Location = TOXL 8164.12.00 06/09
08/26/2009 11:56

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	11:57
02 Std. Sol.	0.289	11:58
03 Room Air	0.000	11:59
04 Std. Sol.	0.291	12:00
05 Room Air	0.000	12:00
06 Std. Sol.	0.292	12:01
07 Room Air	0.000	12:02


08 Sim Temp = 34.0°C

Simul Ser No = 1

Std Sol No = DR00000

County = 08

Oper No. = 107501

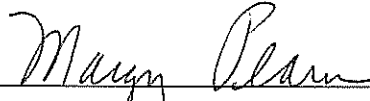

Operator Signature
MARGY PEARSON

Remarks:

Linearity

Form 106-I8000


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-002669
Location = TOXL 8164.12.00 06/09
08/26/2009 14:10

Test	AC	Time
01 Room Air	0.000	14:12
02 *Subject Test	INT*	14:12
03 Room Air	0.000	14:13

*Invalid Test
Interferent Detected

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 11/01/1982
Sub Sex = Female Weight = 150
Test = DUI Cit = NA
Dr. Lic. = ND/DIS821456
Simul Ser No = 1
Std Sol No = DR00000
County = 08 Oper No. = 107501

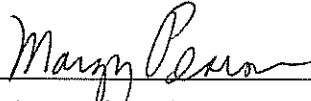

Operator Signature
MARGY PEARSON

Remarks:

ACC tone

Form 106-I8000

*P. 9/2/16
vrm*


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-002669
Location = TOXL 8164.12.00 06/09
08/26/2009 14:00

Test	AC	Time
01 Room Air	0.000	14:01
02 *Subject Test	INT*	14:01
03 Room Air	0.000	14:02

*Invalid Test
Interferent Detected

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 11/01/1982
Sub Sex = Female Weight = 150
Test = DUI Cit = NA
Dr. Lic. = ND/DIS821456
Simul Ser No = 1
Std Sol No = DR00000
County = 08 Oper No. = 107501

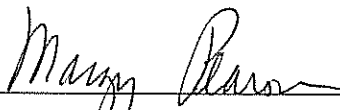

Operator Signature
MARGY PEARSON

Remarks:

Methanol

Form 106-I8000

*P. 10/16/16
JFM*


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-002669
Location = TOXL 8164.12.00 06/09
08/26/2009 14:07

Test	AC	Time
01 Room Air	0.000	14:08
02 *Subject Test	INT*	14:09
03 Room Air	0.000	14:09

*Invalid Test
Interferent Detected

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 11/01/1982
Sub Sex = Female Weight = 150
Test = DUI Cit = NA
Dr. Lic. = ND/DIS821456
Simul Ser No = 1
Std Sol No = DR00000
County = 08 Oper No. = 107501

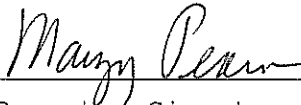

Operator Signature
MARGY PEARSON

Remarks:

1 PA

Form 106-I8000

P. Hoff
pm


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-002669
Location = TOXL 8164.12.00 06/09
08/26/2009 14:19

Test	AC	Time
01 Room Air	AMB*	14:20
02 Room Air	0.000	14:21

*Invalid Test
Check Ambient Conditions

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 11/01/1982
Sub Sex = Female Weight = 150
Test = DUI Cit = NA
Dr. Lic. = ND/DIS821456
Simul Ser No = 1
Std Sol No = DR00000
County = 08 Oper No. = 107501

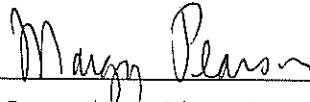

Operator Signature
MARGY PEARSON

Remarks:

Ambient

Form 106-I8000

*P. 12 of 14
PM*


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-002669
Location = TOXL 8164.12.00 06/09
08/26/2009 14:22

Test	AC	Time
01 Room Air	0.000	14:23
02 Invalid Sample	X.XXX	14:23
03 Room Air	0.000	14:24

*Invalid Test - Mouth Alcohol

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 11/01/1982
Sub Sex = Female Weight = 150
Test = DUI Cit = NA
Dr. Lic. = ND/DIS821456
Simul Ser No = 1
Std Sol No = DR00000
County = 08 Oper No. = 107501

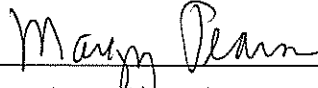

Operator Signature
MARGY PEARSON

Remarks:

Invalid Sample
Form 106-18000

P.13 of 14
JK

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501


CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-002669
Location = TOXL 8164.12.00 06/09
08/26/2009 09:34

Test	AC	Time
01 Diagnostic	OK	09:35
02 Room Air	RFI*	09:35
03 Room Air	0.000	09:36

*Invalid Test
Inhibited - RFI

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 11/01/1982
Sub Sex = Female Weight = 150
Test = DUI Cit = NA
Dr. Lic. = ND/DIS821456
Simul Ser No = 1
Std Sol No = DR00000
County = 08 Oper No. = 107501

I followed the Approved Method and the instructions displayed by the Intoxilyzer in conducting this test.

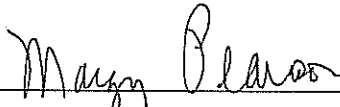

Operator Signature
MARGY PEARSON

Remarks:

RFI

Form 106-I8000

P. 7/24/14
FM

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-002669
Location = TOXL 8164.12.00 06/09
08/26/2009 14:36

Test	AC	Time
01 Room Air	0.000	14:37
02 Std. Gas	0.081	14:37
03 Room Air	0.000	14:38
04 Std. Gas	0.080	14:38
05 Room Air	0.000	14:39
06 Std. Gas	0.081	14:39
07 Room Air	0.000	14:39
08 Std. Gas	0.081	14:40
09 Room Air	0.000	14:40
10 Std. Gas	0.081	14:41
11 Room Air	0.000	14:41
12 Std. Gas	0.080	14:42
13 Room Air	0.000	14:42
14 Std. Gas	0.081	14:42
15 Room Air	0.000	14:43
16 Std. Gas	0.081	14:43
17 Room Air	0.000	14:44
18 Std. Gas	0.081	14:44
19 Room Air	0.000	14:45
20 Std. Gas	0.081	14:45
21 Room Air	0.000	14:46
22 Std. Gas	0.081	14:46
23 Room Air	0.000	14:46
24 Std. Gas	0.082	14:47
25 Room Air	0.000	14:47
26 Std. Gas	0.080	14:48
27 Room Air	0.000	14:48
28 Std. Gas	0.081	14:48
29 Room Air	0.000	14:49
30 Std. Gas	0.081	14:49
31 Room Air	0.000	14:50
32 Std. Gas	0.080	14:50
33 Room Air	0.000	14:51
34 Std. Gas	0.081	14:51
35 Room Air	0.000	14:52
36 Std. Gas	0.081	14:52
37 Room Air	0.000	14:52
38 Std. Gas	0.081	14:53
39 Room Air	0.000	14:53
40 Std. Gas	0.081	14:54
41 Room Air	0.000	14:54
42 Cal Check Stats:		
43 Average	0.0809	
44 Std Dev	0.0005	
45 Rel Std Dev	0.61%	

Stability #1 of 2
MAP

P.15 of 14
RM

Operator Signature

MARGY PEARSON

Remarks:

Lot No 617774

Cyl No 32

Exp Date 5/28/2011

Form 106-I8000

Stability #2 of 2

MAP

P. Kopf/4
15/11