

✓ CML

INTOXILYZER 5000 G and 8000 INSPECTION

Serial No.: 80-003060

Location: TOX L

Check When Done

- A. Is the warm-up time less than 20 minutes? ✓
- B. Is a three-pronged grounded outlet used? ✓
- C. Is the breath tube heated? ✓
- D. Is the diagnostic check complete? ✓
- E. Is the time, date, and year correct (re-set if necessary)? ✓
- F. Print test. (Attach test record.) ✓

---

- G. Low AC. Use <0.03 AC in ACA mode. (Attach test record.) ✓  
Sim SN DR5190 Lot No. — Rep AC 0.021
- H. Linearity Test. Use ≥ 0.25 AC in ACA mode. (Attach test record.) ✓  
Sim SN DR5113 Lot No. — Rep AC 0.297
- I. Interferent check. Use 0.05% acetone plus 0.10 AC ethanol in ABA mode. (Attach test record.) ✓  
Sim SN DR5135 Lot No. ✓ Rep AC 0.10 AC + Acetone
- J. RFI check. Run CMS mode. Key radio on first room air. (Attach test record.) ✓

K. Calibration Check. Use Ethanol Breath Standard cylinder. Do three sets in ACA mode. (Attach test record.)

Lot No. 622252 Cyl No. 28 Exp. Date 5/26/2011

Test 1 <u>0.079</u>	Test 4 <u>0.080</u>	Test 7 <u>0.080</u>
Test 2 <u>0.080</u>	Test 5 <u>0.080</u>	Test 8 <u>0.080</u>
Test 3 <u>0.079</u>	Test 6 <u>0.080</u>	Test 9 <u>0.080</u>

Average 0.080

L. Remarks/Maintenance record: Annual Inspection

27 Aug 2009  
Date

04 Sept 2009  
Date

Mary Ploran  
Field Inspector's Signature

Deb Kashur  
Reviewed by



Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501


CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-003060  
Location = TOXL      8164.12.00 06/09  
08/27/2009      10:33

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	10:33
02 Std. Sol.	0.019	10:39
03 Room Air	0.000	10:40
04 Std. Sol.	0.019	10:40
05 Room Air	0.000	10:41
06 Std. Sol.	0.020	10:42
07 Room Air	0.000	10:42
08 Sim Temp = 34.0°C		

Simul Ser No = DR00000  
Std Sol No = 0000  
County = 08

Oper No. = 107501



Operator Signature  
MARGY PEARSON

Remarks:

LOW 0.021 AC

Form 106-I8000

80-003060  
P. 3 of 9 JAH

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-003060  
Location = TOXL      8164.12.00 06/09  
08/27/2009      11:02

WET CAL CHECK


Test	AC	Time
01 Room Air	0.000	11:03
02 Std. Sol.	0.295	11:05
03 Room Air	0.000	11:05
04 Std. Sol.	0.301	11:06
05 Room Air	0.000	11:07
06 Std. Sol.	0.300	11:07
07 Room Air	0.000	11:08
08 Sim Temp = 34.0°C		

Simul Ser No = DR00000

Std Sol No = 0000

County = 08

Oper No. = 107501

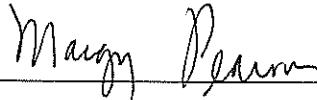

Operator Signature  
MARGY PEARSON

Remarks:

Linearity 0.297 AC

Form 106-I8000

80-003060  
P. 4 of 9 PK


Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-003060  
Location = TOXL      8164.12.00 06/09  
08/27/2009      11:09

Test	AC	Time
01 Room Air	0.000	11:10
02 *Subject Test	INT*	11:10
03 Room Air	0.000	11:11

\*Invalid Test  
Interferent Detected

Sub Name = DISCOVER, THE SPIRIT  
Sub DOB = 11/01/1982  
Sub Sex = Female      Weight = 150  
Test = DUI      Cit = NA  
Dr. Lic. = ND/DIS821456  
Simul Ser No = DR00000  
Std Sol No = 0000  
County = 08      Oper No. = 107501

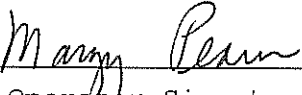

Operator Signature  
MARGY PEARSON

Remarks:

*Interfer*

Form 106-I8000

80-003060  
P. 5 of 9 JH

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-003060  
Location = TOXL      8164.12.00 06/09  
08/27/2009      09:41

Test	AC	Time
01 Diagnostic	OK	09:42
02 Room Air	RFI*	09:42
03 Room Air	0.000	09:42

\*Invalid Test  
Inhibited - RFI

Sub Name = DISCOVER, THE SPIRIT  
Sub DOB = 11/01/1982  
Sub Sex = Female      Weight = 150  
Test = DUI      Cit = NA  
Dr. Lic. = ND/DIS821456  
Lot No = 622252  
Cyl No = 28  
Expiration Date = 05/28/2011  
County = 08      Oper No. = 107501

I followed the Approved Method and the instructions displayed by the Intoxilyzer in conducting this test.



Operator Signature  
MARGY PEARSON

Remarks:

RFI

Form 106-I8000

80-003060  
P. 6 of 9 1/2

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-003060  
Location = TOXL      8164.12.00 06/09  
08/27/2009      09:19

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	09:19
02 Std. Gas	0.079	09:20
03 Room Air	0.000	09:20
04 Std. Gas	0.080	09:20
05 Room Air	0.000	09:21
06 Std. Gas	0.079	09:21
07 Room Air	0.000	09:22

Lot No = 622252

Cyl No = 28

Exp Date = 05/28/2011

County = 08

Oper No. = 107501



Operator Signature  
MARGY PEARSON

Remarks:

ACA

Form 106-I8000

80-003060  
P. 2 of 9 YAK


Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-003060  
Location = TOXL      8164.12.00 06/09  
08/27/2009      09:28

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	09:29
02 Std. Gas	0.080	09:29
03 Room Air	0.000	09:30
04 Std. Gas	0.080	09:30
05 Room Air	0.000	09:31
06 Std. Gas	0.080	09:31
07 Room Air	0.000	09:31

Lot No = 622252  
Cyl No = 28  
Exp Date = 05/28/2011  
County = 08      Oper No. = 107501

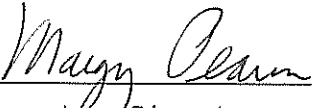  
\_\_\_\_\_  
Operator Signature  
MARGY PEARSON

Remarks:

ACA

Form 106-I8000

80-003060  
P. 8 of 9 AL



Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-003060  
Location = TOXL      8164.12.00 06/09  
08/27/2009      09:23

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	09:23
02 Std. Gas	0.080	09:24
03 Room Air	0.000	09:24
04 Std. Gas	0.080	09:25
05 Room Air	0.000	09:25
06 Std. Gas	0.080	09:26
07 Room Air	0.000	09:26


Lot No = 622252

Cyl No = 28

Exp Date = 05/28/2011

County = 08

Oper No. = 107501

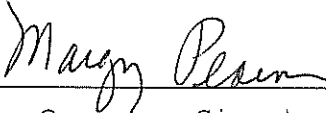

Operator Signature  
MARGY PEARSON

Remarks:

ACA

Form 106-I8000

80-003060  
P. 9 of 9 AL