

✓CML

INTOXILYZER 5000 G and 8000 INSPECTION

Serial No.: 80-004191

Location: TOXL

Check When Done

- A. Is the warm-up time less than 20 minutes? ✓
- B. Is a three-pronged grounded outlet used? ✓
- C. Is the breath tube heated? ✓
- D. Is the diagnostic check complete? ✓
- E. Is the time, date, and year correct (re-set if necessary)? ✓
- F. Print test. (Attach test record.) ✓
- G. Low AC. Use <0.03 AC in ACA mode. (Attach test record.) ✓
Sim SN DR5190 Lot No. — Rep AC 0.021
- H. Linearity Test. Use ≥ 0.25 AC in ACA mode. (Attach test record.) ✓
Sim SN DR 5113 Lot No. — Rep AC 0.297
- I. Interferent check. Use 0.05% acetone plus 0.10 AC ethanol in ABA mode. (Attach test record.) ✓
Sim SN DR5135 Lot No. — Rep AC 0.10 AC + Acetone
- J. RFI check. Run CMS mode. Key radio on first room air. (Attach test record.) ✓
- K. Calibration Check. Use Ethanol Breath Standard cylinder. Do three sets in ACA mode. (Attach test record.)

Lot No. 617774 Cyl No. 34 Exp. Date 5/28/2011

Test 1 <u>0.078</u>	Test 4 <u>0.078</u>	Test 7 <u>0.079</u>
Test 2 <u>0.079</u>	Test 5 <u>0.078</u>	Test 8 <u>0.079</u>
Test 3 <u>0.079</u>	Test 6 <u>0.079</u>	Test 9 <u>0.079</u>

Average 0.079 AC

L. Remarks/Maintenance record: Annual Inspection

1 September 2009
Date

04 Sept 2009
Date

Mangy Peem
Field Inspector's Signature

Deb Kashur
Reviewed by

08/09

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501


CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004191
Location = TOXL 8164.12.00 06/09
09/01/2009 16:29

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	16:30
02 Std. Sol.	0.016	16:30
03 Room Air	0.000	16:31
04 Std. Sol.	0.018	16:32
05 Room Air	0.000	16:32
06 Std. Sol.	0.017	16:33
07 Room Air	0.000	16:34

08 Sim Temp = 34.0°C

Simul Ser No = DR0000
Std Sol No = 0000
County = 08 Oper No. = 107501

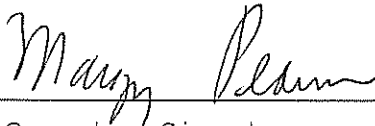

Operator Signature
MARGY PEARSON

Remarks:

Low 0.021 AC

Form 106-I8000

80-004191
P. 3 of 9 FW


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004191
Location = TOXL 8164.12.00 06/09
09/01/2009 16:43

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	16:44
02 Std. Sol.	0.299	16:45
03 Room Air	0.000	16:46
04 Std. Sol.	0.302	16:46
05 Room Air	0.000	16:47
06 Std. Sol.	0.301	16:48
07 Room Air	0.000	16:48
08 Sim Temp = 34.0°C		

Simul Ser No = DR0000
Std Sol No = 0000
County = 08 Oper No. = 107501

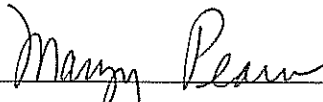

Operator Signature
MARGY PEARSON

Remarks:

Linearity

Form 106-I8000

80-004191
P. 4 of 8 FILE


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004191
Location = TOXL 8164.12.00 06/09
09/01/2009 17:19

Test	AC	Time
01 Room Air	0.000	17:20
02 *Subject Test	INT*	17:21
03 Room Air	0.000	17:21

*Invalid Test
Interferent Detected

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 11/01/1982
Sub Sex = Female Weight = 150
Test = DUI Cit = NA
Dr. Lic. = ND/DIS821456
Simul Ser No = DROOOO
Std Sol No = 0000
County = 08 Oper No. = 107501

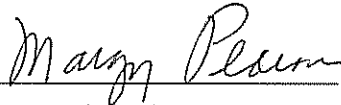

Operator Signature
MARGY PEARSON

Remarks:

Alc + Acetone

Form 106-18000

*80-004191
p. 5 of 9 pr*


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004191
Location = TOXL 8164.12.00 06/09
09/01/2009 15:55

Test	AC	Time
01 Room Air	RFI*	15:56
02 Room Air	0.000	15:56

*Invalid Test
Inhibited - RFI

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 11/01/1982
Sub Sex = Female Weight = 150
Test = DUI Cit = NA
Dr. Lic. = ND/DIS821456
Simul Ser No = DROOOO
Std Sol No = 0000
County = 08 Oper No. = 107501

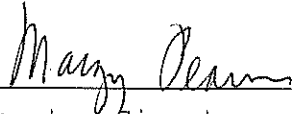

Operator Signature
MARGY PEARSON

Remarks:

RFI

Form 106-18000

80-004191
P. 6 of 9 pm

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004191
Location = TOXL 8164.12.00 06/09
09/01/2009 13:40

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	13:40
02 Std. Gas	0.078	13:41
03 Room Air	0.000	13:41
04 Std. Gas	0.079	13:42
05 Room Air	0.000	13:42
06 Std. Gas	0.079	13:43
07 Room Air	0.000	13:43


Lot No = 617774

Cyl No = 34

Exp Date = 05/28/2011

County = 08

Oper No. = 107501

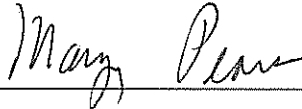

Operator Signature
MARGY PEARSON

Remarks:

ACA

Form 106-I8000

80-004191
p. 7 of 9 *PL*

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004191
Location = TOXL 8164.12.00 06/09
09/01/2009 13:46

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	13:46
02 Std. Gas	0.078	13:47
03 Room Air	0.000	13:47
04 Std. Gas	0.078	13:48
05 Room Air	0.000	13:48
06 Std. Gas	0.079	13:49
07 Room Air	0.000	13:49


Lot No = 617774

Cyl No = 34

Exp Date = 05/28/2011

County = 08

Oper No. = 107501

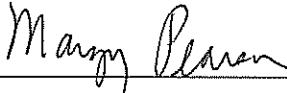

Operator Signature
MARGY PEARSON

Remarks:

ACA

Form 106-I8000

80-004191
P. 8 of 9 PK

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004191
Location = TOXL 8164.12.00 06/09
09/01/2009 13:52

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	13:53
02 Std. Gas	0.079	13:53
03 Room Air	0.000	13:54
04 Std. Gas	0.079	13:54
05 Room Air	0.000	13:55
06 Std. Gas	0.079	13:55
07 Room Air	0.000	13:56


Lot No = 617774

Cyl No = 34

Exp Date = 05/28/2011

County = 08

Oper No. = 107501

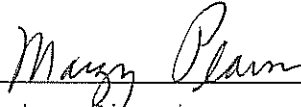

Operator Signature
MARGY PEARSON

Remarks:

ACA

Form 106-I8000

80-004191
8-9-09 JTL