

✓ CML

INTOXILYZER 5000 G and 8000 INSPECTION

Serial No.: 80-004199

Location: TOXL

Check When Done

- A. Is the warm-up time less than 20 minutes? ✓
- B. Is a three-pronged grounded outlet used? ✓
- C. Is the breath tube heated? ✓
- D. Is the diagnostic check complete? ✓
- E. Is the time, date, and year correct (re-set if necessary)? ✓
- F. Print test. (Attach test record.) ✓

- G. Low AC. Use <0.03 AC in ACA mode. (Attach test record.) ✓
Sim SN D125190 Lot No. — Rep AC 0.021
- H. Linearity Test. Use ≥ 0.25 AC in ACA mode. (Attach test record.) ✓
Sim SN DR5113 Lot No. — Rep AC 0.2977%
- I. Interferent check. Use 0.05% acetone plus 0.10 AC ethanol in ABA mode. (Attach test record.) ✓
Sim SN BR 5135 Lot No. — Rep AC 10 AC + Acetone
- J. RFI check. Run CMS mode. Key radio on first room air. (Attach test record.) ✓
- K. Calibration Check. Use Ethanol Breath Standard cylinder. Do three sets in ACA mode. (Attach test record.)

Lot No. 617774 Cyl No. 36 Exp. Date 5/28/2011

Test 1 <u>0.078</u>	Test 4 <u>0.077</u>	Test 7 <u>0.078</u>
Test 2 <u>0.078</u>	Test 5 <u>0.078</u>	Test 8 <u>0.078</u>
Test 3 <u>0.077</u>	Test 6 <u>0.078</u>	Test 9 <u>0.078</u>

Average 0.078

L. Remarks/Maintenance record: Annual Inspection

1 Sept 2009 Date Mary Plam Field Inspector's Signature

04 Sept 2009 Date D. b Kashner Reviewed by

80-004199
P. 1 of 9 FM

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004199
Location = TOXL 8164.12.00 06/09
09/01/2009 14:29

***** Printer Test *****

abcdefghijklmnopqrstuvwxyz1234567890-=_|
ABCDEFGHIJKLMNPOQRSTUVWXYZ!@#\$\$%^&*()_+?

abcdefghijklmnopqrstuvwxyz1234567890-=_|
ABCDEFGHIJKLMNPOQRSTUVWXYZ!@#\$\$%^&*()_+?

Current Instrument Setup

Data Entry Mode: Enabled
Start Test Sequence: DABACABA
Display Prelim Rslt? Yes
Display Third Digit? Yes
Inhib Printer(Y/N)? No
Display Volume? Yes
Disable On Memfull? No

of Print Copies? 1
Select Std (D/W/I)? Dry
Standard Value? 0.080
Standard Lot #? 617774
Standard Cyl #? 36
Standard Expiration? 05/28/2011
Oper No? 107501

Flow Cal. Date: 06/04/2009
Slope 701
Intercept -580043

IR Calibration Date: 06/03/2009
 3um 9um

0th Coef(*100): -44451 -27073
1st Coef(*100): 294251 132474
2nd Coef(*100): -1375 1189
H2O adj(mg/l*10k): 315 545

***** Printer Test End *****



Operator Signature
MARGY PEARSON

Remarks:

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501


CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004199
Location = TOXL 8164.12.00 06/09
09/01/2009 16:45

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	16:45
02 Std. Sol.	0.017	16:46
03 Room Air	0.000	16:47
04 Std. Sol.	0.018	16:47
05 Room Air	0.000	16:48
06 Std. Sol.	0.017	16:49
07 Room Air	0.000	16:49

08 Sim Temp = 34.0°C

Simul Ser No = DR0000
Std Sol No = 000
County = 08 Oper No. = 107501

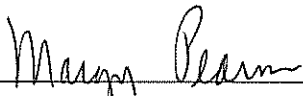

Operator Signature
MARGY PEARSON

Remarks:

Low AC 0.021

Form 106-I8000

80-004199
P. 3-A-9 FH

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004199
Location = TOXL 8164.12.00 06/09
09/01/2009 16:25

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	16:25
02 Std. Sol.	0.299	16:26
03 Room Air	0.000	16:27
04 Std. Sol.	0.299	16:28
05 Room Air	0.000	16:28
06 Std. Sol.	0.298	16:29
07 Room Air	0.000	16:30


08 Sim Temp = 34.0°C

Simul Ser No = DR0000

Std Sol No = 000

County = 08

Oper No. = 107501

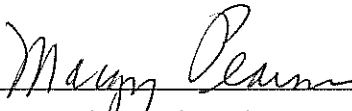

Operator Signature
MARGY PEARSON

Remarks:

Linearity

Form 106-18000

80-004199
P. 4/9 FM


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004199
Location = TOXL 8164.12.00 06/09
09/01/2009 16:53

Test	AC	Time
01 Room Air	0.000	16:53
02 *Subject Test	INT*	16:54
03 Room Air	0.000	16:54

*Invalid Test
Interferent Detected

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 11/01/1982
Sub Sex = Female Weight = 150
Test = DUI Cit = NA
Dr. Lic. = ND/DIS821456
Simul Ser No = DR0000
Std Sol No = 000
County = 08 Oper No. = 107501

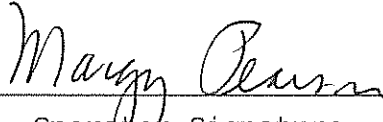

Operator Signature
MARGY PEARSON

Remarks:

Alc + Acetone

Form 106-I8000

80-004199
P. 5 of 9 PL


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004199
Location = TOXL 8164.12.00 06/09
09/01/2009 15:48

Test	AC	Time
01 Room Air	RFI*	15:48
02 Room Air	0.000	15:49

*Invalid Test
Inhibited - RFI

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 11/01/1982
Sub Sex = Female Weight = 150
Test = DUI Cit = NA
Dr. Lic. = ND/DIS821456
Lot No = 617774
Cyl No = 36
Expiration Date = 05/28/2011
County = 08 Oper No. = 107501

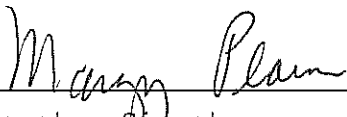

Operator Signature
MARGY PEARSON

Remarks:

RFI

Form 106-I8000

80-004199
P. 6 of 9 FM

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004199
Location = TOXL 8164.12.00 06/09
09/01/2009 13:58

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	13:58
02 Std. Gas	0.078	13:59
03 Room Air	0.000	13:59
04 Std. Gas	0.078	13:59
05 Room Air	0.000	14:00
06 Std. Gas	0.077	14:00
07 Room Air	0.000	14:01


Lot No = 617774

Cyl No = 36

Exp Date = 05/28/2011

County = 08

Oper No. = 107501

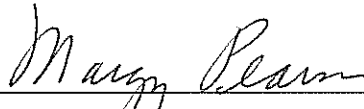

Operator Signature
MARGY PEARSON

Remarks:

ACA

Form 106-I8000

80-004199
P. 7 of 9 RW


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004199
Location = TOXL 8164.12.00 06/09
09/01/2009 14:20

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	14:21
02 Std. Gas	0.077	14:21
03 Room Air	0.000	14:21
04 Std. Gas	0.078	14:22
05 Room Air	0.000	14:22
06 Std. Gas	0.078	14:23
07 Room Air	0.000	14:23

Lot No = 617774
Cyl No = 36
Exp Date = 05/28/2011
County = 08 Oper No. = 107501



Operator Signature
MARGY PEARSON

Remarks:

ACA

Form 106-I8000

80-004199
P. 8 of 9 19h


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004199
Location = TOXL 8164.12.00 06/09
09/01/2009 14:24

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	14:25
02 Std. Gas	0.078	14:25
03 Room Air	0.000	14:26
04 Std. Gas	0.078	14:26
05 Room Air	0.000	14:27
06 Std. Gas	0.078	14:27
07 Room Air	0.000	14:28

Lot No = 617774
Cyl No = 36
Exp Date = 05/28/2011
County = 08 Oper No. = 107501

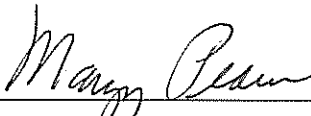
Operator Signature
MARGY PEARSON

Remarks:

ACA

Form 106-I8000

PO-004199
P 9 of 9 FM