

Johnny,

Thank you for sending the Installs back to the Lab so quickly!

However, I'm sending it back to you so you can make a couple corrections.

The cylinder Install looks good, but I need another calibration test (ACA) for the ⁸⁰⁰⁰ Install/Repair Checkout.

Also, I notice you ran the RFI test in ACA mode? I see that because on the RFI form test you sent it it has Dry Cal Check on the top. Please run another RFI test but in CMS mode → Esc Esc Start test and key the mic on the 1st room air.

Then, cross out the date, with one line thru it, on the 8000 ~~Install~~ sheet, put the new ACA date above the old date & initial next to it.

Recap:

1. ACA test
2. RFI test in CMS mode
3. New date on Install form, but you can keep the date on the Cyl. Install form since that one is OK.

Thanks Johnny,
Deb Kashur
328-6159

80-009188
P. 1/16 PM

✓ Database
cm

**INTOXILYZER® 5000 G AND 8000
INSTALLATION AND REPAIR CHECKOUT**

Serial No.: 80-004188 Instrument Location: Bowman Co SO

Reason for Install/Repair: _____ Install After Lab Repair/Inspection
X Other (specify) NEW Instrument

Check When Done:

- X A. Surge Protector Installed/Properly Grounded.
- N/A B. Telephone Line Connected to Intoxilyzer® KB.
- X C. Breath Tube Heated.
- X D. Review/Enter Preliminary Data Entry (Level 2, Function E).
- X E. Ethanol Breath Standard: _____ Connected X Installed.
(Include Form 105-G if New Cylinder Installed.)
- X F. Review/Enter Gas Setup (Level 1, Function S). (Intoxilyzer® 8000 Only.)
- X G. Print Test (Level 1, Function P). (Attach Test Record.)
- X H. ACA Test (Level 1, Function C). (Attach Test Record.)
- X I. RFI Test (Use CMS Mode-First Room Air). (Attach Test Record.)
- J. Repair and/or Maintenance Performed: _____

9-30-09 10-01-09 ^{map}
Date
10-13-09
Date

[Signature]
Field Inspector's Signature
Maple
Reviewed By

80-004188
P. 2 of 6

Intoxilyzer Test Record and Checklist
 NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
 North Dakota Model 8000 SN 80-004188
 Location = BOWM 8164.12.00 06/09
~~09/30/2009~~ 09:07
 10-01-2009 MAB

***** Printer Test *****

abcdefghijklmnopqrstuvwxy1234567890-=_|
 ABCDEFGHIJKLMNOPQRSTUVWXYZ!@#\$\$%^&*()_+?

abcdefghijklmnopqrstuvwxy1234567890-=_|
 ABCDEFGHIJKLMNOPQRSTUVWXYZ!@#\$\$%^&*()_+?

Current Instrument Setup


Data Entry Mode: Enabled
 Start Test Sequence: DABACABA
 Display Prelim Rslt? Yes
 Display Third Digit? Yes
 Inhib Printer(Y/N)? No
 Display Volume? No
 Disable On Memfull? Yes
 # of Print Copies? 1
 Select Std (D/W/I)? Dry
 Standard Value? 0.080
 Standard Lot #? 659358
 Standard Cyl #? 51
 Standard Expiration? 09/03/2011
 Oper No? 020814

Flow Cal. Date: 06/04/2009
 Slope 689
 Intercept -630047

IR Calibration Date: 06/03/2009
 3um 9um

	3um	9um
0th Coef(*100):	-23767	-23243
1st Coef(*100):	280491	137567
2nd Coef(*100):	1648	933
H2O adj(mg/l*10k):	733	647

***** Printer Test End *****

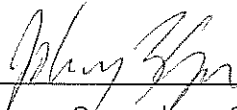

 Operator Signature
 JOHNNY Z. LAWSON

Remarks: Print Test


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004188
Location = BOWM 8164.12.00 06/09
10/01/2009 14:17

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	14:18
02 Std. Gas	0.083	14:18
03 Room Air	0.000	14:19
04 Std. Gas	0.083	14:19
05 Room Air	0.000	14:20
06 Std. Gas	0.084	14:20
07 Room Air	0.000	14:21

Lot No = 659358
Cyl No = 51
Exp Date = 09/03/2011
County = 06 Oper No. = 020814

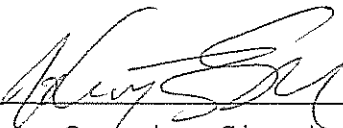
Operator Signature
JOHNNY Z. LAWSON

Remarks:

ACA TEST
Form 106-I8000

80-004188
P. 4-16-10

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004188
Location = BOWM 8164.12.00 06/09
10/01/2009 14:23


Test	AC	Time
01 Diagnostic	OK	14:25
02 Room Air	RFI*	14:25
03 Room Air	0.000	14:26

*Invalid Test
Inhibited - RFI

Sub Name = TEST, SUBJECT NA
Sub DOB = 08/08/1990
Sub Sex = Male Weight = 200
Test = OTH Cit = 00000000
Dr. Lic. = ND/090909
Lot No = 659358
Cyl No = 51

Expiration Date = 09/03/2011
County = 06 Oper No. = 020814

I followed the Approved Method and the instructions displayed by the Intoxilyzer in conducting this test.


Operator Signature
JOHNNY Z. LAWSON

Remarks:

RFI TEST

Form 106-I8000

80-004188
P. Sofle fm

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501


CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004188
Location = BOWM 8164.12.00 06/09
09/30/2009 09:17

DRY CAL CHECK

Test	AC	Time
01 Room Air	RFI*	09:17
02 Room Air	RFI*	09:17

*Invalid Test
Inhibited - RFI

Lot No = 659358
Cyl No = 51
Exp Date = 09/03/2011
County = 06 Oper No. = 020814

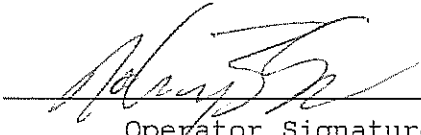

Operator Signature
JOHNNY Z. LAWSON

Remarks:

Form 106-I8000

80-004188
P. 6 of 15