

Database ✓

**INTOXILYZER® 5000 G AND 8000
INSTALLATION AND REPAIR CHECKOUT**

Serial No.: 80-004196 Instrument Location: Williston LEC ^{PD}

Reason for Install/Repair: _____ Install After Lab Repair/Inspection
X Other (specify) New install

Check When Done:

- ✓ A. Surge Protector Installed/Properly Grounded.
- ✓ B. Telephone Line Connected to Intoxilyzer® KB.
- ✓ C. Breath Tube Heated.
- ✓ D. Review/Enter Preliminary Data Entry (Level 2, Function E).
- ✓ E. Ethanol Breath Standard: ✓ Connected ✓ Installed.
(Include Form 105-G if New Cylinder Installed.)
- ✓ F. Review/Enter Gas Setup (Level 1, Function S). (Intoxilyzer® 8000 Only.)
- ✓ G. Print Test (Level 1, Function P). (Attach Test Record.)
- ✓ H. ACA Test (Level 1, Function C). (Attach Test Record.)
- ✓ I. RFI Test (Use CMS Mode-First Room Air). (Attach Test Record.)
- J. Repair and/or Maintenance Performed: _____

10/01/2009
Date
06 Oct 09
Date

Darryl W. Aberle
Field Inspector's Signature
Deb Keshner
Reviewed By

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004196
Location = WILL 8164.12.00 06/09
10/01/2009 17:37

***** Printer Test *****

abcdefghijklmnopqrstuvwxyz1234567890-=_|
ABCDEFGHIJKLMNPOQRSTUVWXYZ!@#\$\$%^&*()_+?

abcdefghijklmnopqrstuvwxyz1234567890-=_|
ABCDEFGHIJKLMNPOQRSTUVWXYZ!@#\$\$%^&*()_+?

Current Instrument Setup


Data Entry Mode: Enabled
Start Test Sequence: DABACABA
Display Prelim Rslt? Yes
Display Third Digit? Yes
Inhib Printer(Y/N)? No
Display Volume? No
Disable On Memfull? Yes
of Print Copies? 1
Select Std (D/W/I)? Dry
Standard Value? 0.080
Standard Lot #? 659357
Standard Cyl #? 45
Standard Expiration? 09/03/2011
Oper No? 010101

Flow Cal. Date: 06/04/2009
Slope 727
Intercept -600779

IR Calibration Date: 06/03/2009
 3um 9um

0th Coef(*100):	-16078	-19878
1st Coef(*100):	252494	128485
2nd Coef(*100):	3196	1141
H2O adj(mg/l*10k):	686	431

***** Printer Test End *****

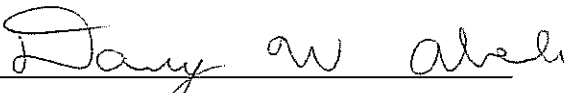

Operator Signature
DARCY W. ABERLE

Remarks:

Printer Test

Williston

Form 106-I8000

80-004196
P. 2 of 4


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004196
Location = WILL 8164.12.00 06/09
10/01/2009 17:38

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	17:38
02 Std. Gas	0.079	17:39
03 Room Air	0.000	17:39
04 Std. Gas	0.081	17:39
05 Room Air	0.000	17:40
06 Std. Gas	0.081	17:40
07 Room Air	0.000	17:41

Lot No = 659357
Cyl No = 45
Exp Date = 09/03/2011
County = 53 Oper No. = 010101

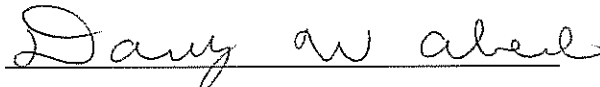
Operator Signature
DARCY W. ABERLE

Remarks:

ACA
Williston

Form 106-18000

80-004196
P.3.4 PM

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501


CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004196
Location = WILL 8164.12.00 06/09
10/01/2009 17:42

Test	AC	Time
01 Diagnostic	OK	17:45
02 Room Air	RFI*	17:45
03 Room Air	0.000	17:46

*Invalid Test
Inhibited - RFI

Sub Name = TEST, SUBJECT N
Sub DOB = 10/01/1980
Sub Sex = Male Weight = 200
Test = OTH Cit = 010101
Dr. Lic. = ND/DIS804455
Lot No = 659357
Cyl No = 45
Expiration Date = 09/03/2011
County = 53 Oper No. = 010101

I followed the Approved Method and the instructions displayed by the Intoxilyzer in conducting this test.

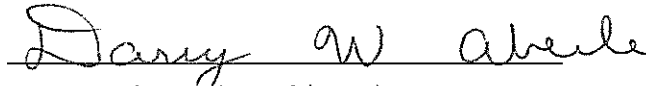

Operator Signature
DARCY W. ABERLE

Remarks:

RFI Test

Williston

Form 106-18000

80-604196
p. 4 of 4 km