

✓ Database will

INTOXILYZER® 5000 G AND 8000
INSTALLATION AND REPAIR CHECKOUT

Serial No.: 80-003071 Instrument Location: Belmont Police Dept

Reason for Install/Repair: _____ Install After Lab Repair/Inspection
X Other (specify) New Install

Check When Done:

- X A. Surge Protector Installed/Properly Grounded.
- X B. Telephone Line Connected to Intoxilyzer® KB.
- X C. Breath Tube Heated.
- X D. Review/Enter Preliminary Data Entry (Level 2, Function E).
- X E. Ethanol Breath Standard: X Connected X Installed.
(Include Form 105-G if New Cylinder Installed.)
- X F. Review/Enter Gas Setup (Level 1, Function S). (Intoxilyzer® 8000 Only.)
- X G. Print Test (Level 1, Function P). (Attach Test Record.)
- X H. ACA Test (Level 1, Function C). (Attach Test Record.)
- X I. RFI Test (Use CMS Mode-First Room Air). (Attach Test Record.)
- J. Repair and/or Maintenance Performed: New Install

10-9-2009
Date
16 Oct 09
Date

[Signature]
Field Inspector's Signature
Deb Kashner
Reviewed By

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-003071
Location = BELC 8164.12.00 06/09
10/09/2009 21:27

***** Printer Test *****

abcdefghijklmnopqrstuvwxyz1234567890-=
ABCDEFGHIJKLMNPOQRSTUVWXYZ!@#\$\$%^&*()_+?

abcdefghijklmnopqrstuvwxyz1234567890-=
ABCDEFGHIJKLMNPOQRSTUVWXYZ!@#\$\$%^&*()_+?

Current Instrument Setup

Data Entry Mode: Enabled
Start Test Sequence: DABACABA
Display Prelim Rslt? Yes
Display Third Digit? Yes
Inhib Printer(Y/N)? No
Display Volume? No
Disable On Memfull? Yes
of Print Copies? 1
Select Std (D/W/I)? Dry
Standard Value? 0.080
Standard Lot #? 659357
Standard Cyl #? 38
Standard Expiration? 09/03/2011
Oper No? 119802

Flow Cal. Date: 04/08/2008
Slope 696
Intercept -591525

IR Calibration Date: 06/05/2008
3um 9um

0th Coef(*100): -7866 -4732
1st Coef(*100): 276944 136843
2nd Coef(*100): 1691 912
H2O adj(mg/l*10k): 337 201

***** Printer Test End *****



Operator Signature
BRIAN F. CRISSLER

Remarks:

New Fastal


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-003071
Location = BELC 8164.12.00 06/09
10/09/2009 21:50

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	21:51
02 Std. Gas	0.080	21:51
03 Room Air	0.000	21:52
04 Std. Gas	0.081	21:52
05 Room Air	0.000	21:53
06 Std. Gas	0.081	21:53
07 Room Air	0.000	21:54

Lot No = 659357
Cyl No = 38
Exp Date = 09/03/2011
County = 63 Oper No. = 119802


Operator Signature
BRIAN F. CRISSLER

Remarks: *New Install*

Form 106-I8000

80-003071
P. 3 of 4 *FL*


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-003071
Location = BELC 8164.12.00 06/09
10/09/2009 21:44

Test	AC	Time
01 Diagnostic	OK	21:46
02 Room Air	0.000	21:47
03 *Subject Test	RFI*	21:47
04 Room Air	0.000	21:47

*Invalid Test
Inhibited - RFI

I followed the Approved Method and the instructions displayed by the Intoxilyzer in conducting this test.

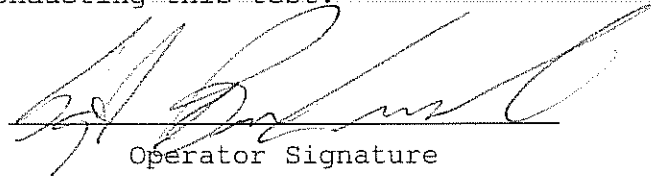

Operator Signature
BRIAN F. CRISSLER

Remarks:

New Install

Form 106-I8000

80-003071
p. 4 of 4 FR