

✓
Database
one

**INTOXILYZER® 5000 G AND 8000
INSTALLATION AND REPAIR CHECKOUT**

Serial No.: 80-003060 Instrument Location: UNDP

Reason for Install/Repair: Install After Lab Repair/Inspection
 Other (specify) _____

Check When Done:

- A. Surge Protector Installed/Properly Grounded.
- B. Telephone Line Connected to Intoxilyzer® KB.
- C. Breath Tube Heated.
- D. Review/Enter Preliminary Data Entry (Level 2, Function E).
- E. Ethanol Breath Standard: Connected Installed.
(Include Form 105-G if New Cylinder Installed.)
- F. Review/Enter Gas Setup (Level 1, Function S). (Intoxilyzer® 8000 Only.)
- G. Print Test (Level 1, Function P). (Attach Test Record.)
- H. ACA Test (Level 1, Function C). (Attach Test Record.)
- I. RFI Test (Use CMS Mode-First Room Air). (Attach Test Record.)
- J. Repair and/or Maintenance Performed: _____

10-14-2009
Date
16 Oct 09
Date

L. Tracy L. Mardinger #258
Field Inspector's Signature
Deb Kasher
Reviewed By

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-003060
Location = UNDP 8164.12.00 06/09
10/14/2009 01:03

***** Printer Test *****

abcdefghijklmnopqrstuvwxyz1234567890-=
ABCDEFGHIJKLMNPOQRSTUVWXYZ!@#%&^*()_+?

abcdefghijklmnopqrstuvwxyz1234567890-=
ABCDEFGHIJKLMNPOQRSTUVWXYZ!@#%&^*()_+?

Current Instrument Setup


Data Entry Mode: Enabled
Start Test Sequence: DABACABA
Display Prelim Rslt? Yes
Display Third Digit? Yes
Inhib Printer(Y/N)? No
Display Volume? No
Disable On Memfull? Yes
of Print Copies? 1
Select Std (D/W/I)? Dry
Standard Value? 0.080
Standard Lot #? 659358
Standard Cyl #? 31
Standard Expiration? 09/03/2011
Oper No? 119409

Flow Cal. Date: 04/08/2008
Slope 702
Intercept -647610

IR Calibration Date: 06/04/2008
 3um 9um

0th Coef(*100): -5190 -7737
1st Coef(*100): 270580 138779
2nd Coef(*100): 2499 956
H2O adj(mg/l*10k): 377 188

***** Printer Test End *****

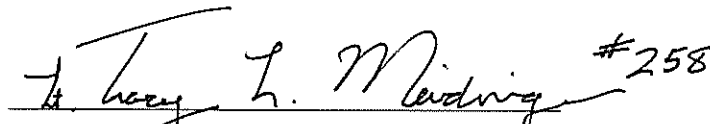 #258

Operator Signature
TRACY L. MEIDINGER

Remarks:

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-003060
Location = UNDP 8164.12.00 06/09
10/14/2009 01:04

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	01:05
02 Std. Gas	0.081	01:05
03 Room Air	0.000	01:06
04 Std. Gas	0.082	01:06
05 Room Air	0.000	01:07
06 Std. Gas	0.082	01:07
07 Room Air	0.000	01:07

Lot No = 659358

Cyl No = 31

Exp Date = 09/03/2011

County = 18

Oper No. = 119409

Tracy L. Meidinger #258
Operator Signature
TRACY L. MEIDINGER

Remarks:

Form 106-I8000

80-003060
p. 3 of 4 FM

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501


CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-003060
Location = UNDP 8164.12.00 06/09
10/14/2009 01:42

Test	AC	Time
01 Diagnostic	OK	01:44
02 Room Air	RFI*	01:44
03 Room Air	0.000	01:45

*Invalid Test
Inhibited - RFI

Sub Name = TEST, SUBJECT GIN
Sub DOB = 09/12/1967
Sub Sex = Male Weight = 190
Test = OTH Cit = 4444444
Dr. Lic. = ND/NA
Lot No = 659358
Cyl No = 31
Expiration Date = 09/03/2011
County = 18 Oper No. = 119409

I followed the Approved Method and the instructions displayed by the Intoxilyzer in conducting this test.

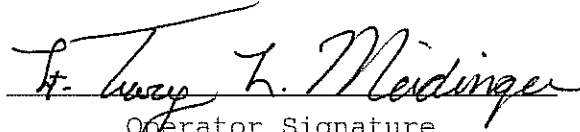

Operator Signature
TRACY L. MEIDINGER

Remarks:

Form 106-I8000

80-003060
P. 4 of 4 JAM