

INTOXILYZER® 8000 INSTALLATION AND REPAIR CHECKOUT

Serial No.: 80-002666 Instrument Location: Williston LEC

Reason for Install/Repair: _____ Install After Lab Repair/Inspection
X Other (Specify) New

Check When Done:

- X A. Surge Protector Installed/Properly Grounded.
- X B. Telephone Line Connected to Intoxilyzer®.
- X C. Breath Tube Heated.
- X D. Review/Enter Preliminary Data Entry (Level 2, Function E).
- X E. Review/Enter Gas Setup (Level 1, Function S).
- X F. Print Test (Level 1, Function P). (Attach Test Record.)
- X G. ACA Test (Level 1, Function C). (Attach Test Record.)
- X H. RFI Test (Use CMS Mode-First Room Air). (Attach Test Record.)
- I. Repair and/or Maintenance Performed: None, New Install

11-18-09
Date

Darryl W. Abene
Field Inspector's Signature

24 NOV 2009
Date

Deb Kashur
Reviewed By

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-002666
Location = WILL 8164.13.00 06/09
11/18/2009 16:18

***** Printer Test *****

abcdefghijklmnopqrstuvwxyz1234567890-=_|
ABCDEFGHIJKLMNPOQRSTUVWXYZ!@#%&^*()_+?

abcdefghijklmnopqrstuvwxyz1234567890-=_|
ABCDEFGHIJKLMNPOQRSTUVWXYZ!@#%&^*()_+?

Current Instrument Setup


Data Entry Mode: Enabled
Start Test Sequence: DABACABA
Display Prelim Rslt? Yes
Display Third Digit? Yes
Inhib Printer(Y/N)? No
Display Volume? No
Disable On Memfull? Yes
of Print Copies? 1
Select Std (D/W/I)? Dry
Standard Value? 0.105
Standard Lot #? 659357
Standard Cyl #? 45
Standard Expiration? 09/03/2011
Oper No? 010101

Flow Cal. Date: 09/10/2007
Slope 677
Intercept -603196

IR Calibration Date: 09/10/2007
 3um 9um

0th Coef(*100): -14148 -8312
1st Coef(*100): 274152 135499
2nd Coef(*100): 1387 855
H2O adj(mg/l*10k): 430 237

***** Printer Test End *****

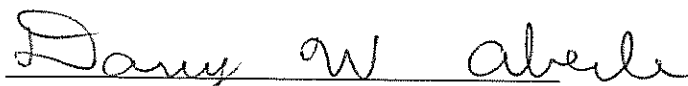

Operator Signature
DARCY W. ABERLE

Remarks:

Print Test


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-002666
Location = WILL 8164.13.00 06/09
11/18/2009 16:18

DRY CAL CHECK

| Test | AC | Time |
|-------------|-------|-------|
| 01 Room Air | 0.000 | 16:19 |
| 02 Std. Gas | 0.081 | 16:19 |
| 03 Room Air | 0.000 | 16:20 |
| 04 Std. Gas | 0.081 | 16:20 |
| 05 Room Air | 0.000 | 16:20 |
| 06 Std. Gas | 0.082 | 16:21 |
| 07 Room Air | 0.000 | 16:21 |

Lot No = 659357
Cyl No = 45
Exp Date = 09/03/2011
County = 53 Oper No. = 010101

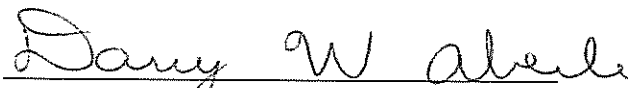

Operator Signature
DARCY W. ABERLE

Remarks:

ACA Test

Form 106-I8000


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-002666
Location = WILL 8164.13.00 06/09
11/18/2009 16:24

| Test | AC | Time |
|-------------|-------|-------|
| 01 Room Air | RFI* | 16:25 |
| 02 Room Air | 0.000 | 16:26 |

*Invalid Test
Inhibited - RFI

Sub Name = TEST, TEST T
Sub DOB = 05/05/1990
Sub Sex = Male Weight = 150
Test = OTH Cit = 000
Dr. Lic. = ND/NA
Lot No = 659357
Cyl No = 45
Expiration Date = 09/03/2011
County = 53 Oper No. = 010101

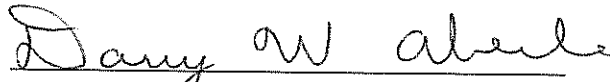

Operator Signature
DARCY W. ABERLE

Remarks:

RFI Test

Form 106-I8000