

INTOXILYZER® 8000 INSPECTION

Serial No.: 80-004936

Location: TOXL

- Check When Done
- A. Is the breath tube heated? ✓
  - B. Is the diagnostic check complete? ✓
  - C. Is the time, date, and year correct (re-set if necessary)? ✓
  - D. Print test. (Attach test record.) ✓
  - E. Low AC. Use  $\leq 0.03$  AC in ACA mode. (Attach test record.)  
Sim SN MP3060 Lot No. 201507F Rep AC 0.020 ✓
  - F. High AC. Use  $\geq 0.25$  AC in ACA mode. (Attach test record.)  
Sim SN MP3070 Lot No. 16060 Rep AC 0.300 ✓
  - G. Interferent check. Use 0.05% acetone in ethanol solution in ABA mode.  
(Attach test record.)  
Sim SN DR5135 Lot No. ICS4 Solution 0.10AC + 0.05% ACETONE ✓
  - H. RFI check. Run CMS mode. Key radio during CMS test.  
(Attach test record.) ✓
  - I. Calibration Check. Use Ethanol Breath Standard cylinder. Do three sets in ACA mode.  
(Attach test record.) 0.080 AC

Lot No. 17316080A2 Cyl No. 026 Exp. Date 8/5/2018  
Test 1 0.080 Test 4 0.080 Test 7 0.081  
Test 2 0.080 Test 5 0.081 Test 8 0.080  
Test 3 0.080 Test 6 0.080 Test 9 0.081  
Average 0.080

J. Tank pressure, Level 3 Func D Sub G P: 543 psi Reg 550 psi

K. Remarks/Maintenance record: ANNUAL INSPECTION, REPLACED  
SIMULATOR RETURN O-RING.

4.4.2017  
Date

Charles E. Egan  
Field Inspector's Signature

04 Apr 17  
Date

Deb Shanaver  
Reviewed by



Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501


CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-004936  
Location = TOXL      8164.14.00 09/16  
04/04/2017      10:53

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	10:53
02 Std. Sol.	0.018	10:54
03 Room Air	0.000	10:55
04 Std. Sol.	0.018	10:55
05 Room Air	0.000	10:56
06 Std. Sol.	0.018	10:57
07 Room Air	0.000	10:57

08 Sim Temp = 34.0°C

Simul Ser No = MP3060  
Std Sol No = 201507F  
County = 08      Oper No. = 666666

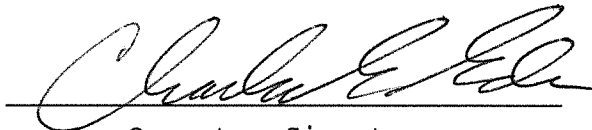

Operator Signature  
CHARLES EDER

Remarks:

*Low AC*  
*0.020 AC*

Form 106-I8000

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer            Alcohol Analyzer  
North Dakota Model 8000            SN 80-004936  
Location = TOXL                    8164.14.00 09/16  
04/04/2017                            11:04

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	11:04
02 Std. Sol.	0.298	11:05
03 Room Air	0.000	11:05
04 Std. Sol.	0.300	11:06
05 Room Air	0.000	11:07
06 Std. Sol.	0.300	11:07
07 Room Air	0.000	11:08

08 Sim Temp = 34.0°C

Simul Ser No = MP3070  
Std Sol No = 16060  
County = 08                            Oper No. = 666666



Operator Signature  
CHARLES EDER

Remarks:    *HIGH AC*  
                  *0.300 AC*

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer            Alcohol Analyzer  
North Dakota Model 8000        SN 80-004936  
Location = TOXL                8164.14.00 09/16  
04/04/2017                        11:09

Test	AC	Time
01 Room Air	0.000	11:10
02 *Subject Test	INT*	11:10
03 Room Air	0.000	11:11

\*Invalid Test  
Interferent Detected

Sub Name = DISCOVER, THE SPIRIT  
Sub DOB = 01/01/1982  
Sub Sex = Female                Weight = 150  
Test = DUI                        Cit = INTERFERENT  
Dr. Lic. = ND/DIS821456  
Lot No = 17316080A2  
Cyl No = 26  
Expiration Date = 08/05/2018  
County = 08                        Oper No. = 666666



Operator Signature  
CHARLES EDER

Remarks:

*Interferent*

*0.10 AC + 0.05% ACETONE*

Form 106-I8000



Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer            Alcohol Analyzer  
North Dakota Model 8000        SN 80-004936  
Location = TOXL                    8164.14.00 09/16  
04/04/2017                            11:15

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	11:16
02 Std. Gas	0.080	11:16
03 Room Air	0.000	11:17
04 Std. Gas	0.080	11:17
05 Room Air	0.000	11:17
06 Std. Gas	0.080	11:18
07 Room Air	0.000	11:18

Lot No = 17316080A2  
Cyl No = 26  
Exp Date = 08/05/2018  
County = 08                            Oper No. = 666666



Operator Signature  
CHARLES EDER

Remarks: CALIBRATION CHECK  
0.080 AC

Form 106-I8000

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-004936  
Location = TOXL      8164.14.00 09/16  
04/04/2017      11:19

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	11:19
02 Std. Gas	0.080	11:20
03 Room Air	0.000	11:20
04 Std. Gas	0.081	11:20
05 Room Air	0.000	11:21
06 Std. Gas	0.080	11:21
07 Room Air	0.000	11:22

Lot No = 17316080A2  
Cyl No = 26  
Exp Date = 08/05/2018  
County = 08      Oper No. = 666666



Operator Signature  
CHARLES EDER

Remarks: CALIBRATION CHECK  
0.080 AC

Form 106-I8000




Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-004936  
Location = TOXL      8164.14.00 09/16  
04/04/2017      11:22

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	11:23
02 Std. Gas	0.081	11:23
03 Room Air	0.000	11:24
04 Std. Gas	0.080	11:24
05 Room Air	0.000	11:25
06 Std. Gas	0.081	11:25
07 Room Air	0.000	11:26

Lot No = 17316080A2  
Cyl No = 26  
Exp Date = 08/05/2018  
County = 08      Oper No. = 666666



Operator Signature  
CHARLES EDER

Remarks: CALIBRATION CHECK  
0.080 AC