

INTOXILYZER® 8000 INSPECTION

Serial No.: 80-004204

Location: TOXL

- | | Check When Done |
|--|-----------------|
| A. Is the breath tube heated? | <u>✓</u> |
| B. Is the diagnostic check complete? | <u>✓</u> |
| C. Is the time, date, and year correct (re-set if necessary)? | <u>✓</u> |
| D. Print test. (Attach test record.) | <u>✓</u> |
| E. Low AC. Use ≤ 0.03 AC in ACA mode. (Attach test record.)
Sim SN <u>MP3003</u> Lot No. <u>201610A</u> Rep AC <u>0.020</u> | <u>✓</u> |
| F. High AC. Use ≥ 0.25 AC in ACA mode. (Attach test record.)
Sim SN <u>MP3063</u> Lot No. <u>16102</u> Rep AC <u>0.400</u> | <u>✓</u> |
| G. Interferent check. Use 0.05% acetone in ethanol solution in ABA mode.
(Attach test record.)
Sim SN <u>DR5135</u> Lot No. <u>ICS4</u> Solution <u>0.10AC + 0.05%
ACETONE</u> | <u>✓</u> |
| H. RFI check. Run CMS mode. Key radio during CMS test.
(Attach test record.) | <u>✓</u> |
| I. Calibration Check. Use Ethanol Breath Standard cylinder. Do three sets in ACA mode.
(Attach test record.) <u>0.080</u> AC | |

Lot No. <u>17316080A2</u>	Cyl No. <u>28</u>	Exp. Date <u>8/5/2018</u>
Test 1 <u>0.079</u>	Test 4 <u>0.080</u>	Test 7 <u>0.080</u>
Test 2 <u>0.080</u>	Test 5 <u>0.080</u>	Test 8 <u>0.080</u>
Test 3 <u>0.080</u>	Test 6 <u>0.079</u>	Test 9 <u>0.080</u>
Average <u>0.080</u>		

J. Tank pressure, Level 3 Func D Sub G P: 786 psi Reg 800 psi

K. Remarks/Maintenance record: REPLACED BREATH HOSE AND SIMULATOR
RETURN O-RINGS, ANNUAL INSPECTION

4-19-17
Date

Charles E. Egan
Field Inspector's Signature

19 Apr '17
Date

Deb Shanover
Reviewed by

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004204
Location = TOXL 8164.14.00 09/16
04/19/2017 14:39

***** Printer Test *****

abcdefghijklmnopqrstuvwxy1234567890-|=|
ABCDEFGHIJKLMNopQRSTUVWXYZ!@#%&*()_+?

abcdefghijklmnopqrstuvwxy1234567890-|=|
ABCDEFGHIJKLMNopQRSTUVWXYZ!@#%&*()_+?

Current Instrument Setup

Data Entry Mode: Enabled
Start Test Sequence: DABACABA
Display Prelim Rslt? Yes
Display Third Digit? Yes
Inhib Printer(Y/N)? No
Display Volume? No
Disable On Memfull? Yes
of Print Copies? 1
Select Std (D/W/I)? Dry
Standard Value? 0.080
Standard Lot #? 17316080A2
Standard Cyl #? 28
Standard Expiration? 08/05/2018
Oper No? 666666

Flow Cal. Date: 07/13/2010
Slope 750
Intercept -815939

IR Calibration Date: 05/12/2015
 3um 9um

0th Coef(*100): -23668 -20605
1st Coef(*100): 259586 131325
2nd Coef(*100): 3023 1492
H2O adj (mg/l*10k): 445 418

***** Printer Test End *****



Operator Signature
CHARLES EDER

Remarks:

PRINT TEST

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004204
Location = TOXL 8164.14.00 09/16
04/19/2017 14:43

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	14:44
02 Std. Sol.	0.019	14:44
03 Room Air	0.000	14:45
04 Std. Sol.	0.019	14:46
05 Room Air	0.000	14:46
06 Std. Sol.	0.019	14:47
07 Room Air	0.000	14:47

08 Sim Temp = 34.0°C

Simul Ser No = MP3003
Std Sol No = 201610A
County = 08 Oper No. = 666666



Operator Signature
CHARLES EDER

Remarks:

Low AC
0.020 AC

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004204
Location = TOXL 8164.14.00 09/16
04/19/2017 14:49

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	14:50
02 Std. Sol.	0.398	14:50
03 Room Air	0.000	14:51
04 Std. Sol.	0.399	14:52
05 Room Air	0.000	14:52
06 Std. Sol.	0.397	14:53
07 Room Air	0.000	14:54

08 Sim Temp = 34.0°C

Simul Ser No = MP3063
Std Sol No = 16102
County = 08 Oper No. = 666666



Operator Signature
CHARLES EDER

Remarks:

HIGH AC
0.400 AC


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004204
Location = TOXL 8164.14.00 09/16
04/19/2017 14:57

Test	AC	Time
01 Room Air	0.000	14:58
02 *Subject Test	INT*	14:59
03 Room Air	0.000	14:59

*Invalid Test
Interferent Detected

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 01/01/1982
Sub Sex = Female Weight = 150
Test = DUI Cit = INTERFERENT
Dr. Lic. = ND/DIS821456
Lot No = 17316080A2
Cyl No = 28
Expiration Date = 08/05/2018
County = 08 Oper No. = 666666

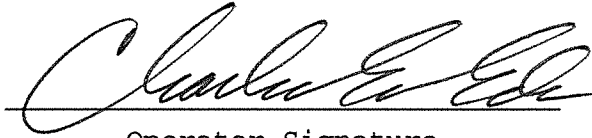

Operator Signature
CHARLES EDER

Remarks: *Interferent*
0.10AC + 0.05% ACETONE

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004204
Location = TOXL 8164.14.00 09/16
04/19/2017 15:00

Test	AC	Time
01 Diagnostic	OK	15:01
02 Room Air	0.000	15:01
03 *Subject Test	RFI*	15:02
04 Room Air	0.000	15:02

*Invalid Test
 Inhibited - RFI

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 01/01/1982
Sub Sex = Female Weight = 150
Test = DUI Cit = RFI CHECK
Dr. Lic. = ND/DIS821456
Lot No = 17316080A2
Cyl No = 28
Expiration Date = 08/05/2018
County = 08 Oper No. = 666666

I followed the Approved Method and the instructions displayed by the Intoxilyzer in conducting this test.



Operator Signature
CHARLES EDER

Remarks:

RFI CHECK


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004204
Location = TOXL 8164.14.00 09/16
04/19/2017 15:06

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	15:06
02 Std. Gas	0.079	15:06
03 Room Air	0.000	15:07
04 Std. Gas	0.080	15:07
05 Room Air	0.000	15:08
06 Std. Gas	0.080	15:08
07 Room Air	0.000	15:09

Lot No = 17316080A2
Cyl No = 28
Exp Date = 08/05/2018
County = 08 Oper No. = 666666

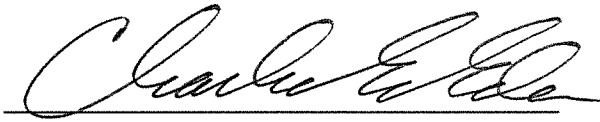

Operator Signature
CHARLES EDER

Remarks: CALIBRATION CHECK
0.080 AC

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004204
Location = TOXL 8164.14.00 09/16
04/19/2017 15:09

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	15:10
02 Std. Gas	0.080	15:10
03 Room Air	0.000	15:10
04 Std. Gas	0.080	15:11
05 Room Air	0.000	15:11
06 Std. Gas	0.079	15:12
07 Room Air	0.000	15:12

Lot No = 17316080A2
Cyl No = 28
Exp Date = 08/05/2018
County = 08 Oper No. = 666666



Operator Signature
CHARLES EDER

Remarks: CALIBRATION CHECK
0.080 AC


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004204
Location = TOXL 8164.14.00 09/16
04/19/2017 15:13

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	15:13
02 Std. Gas	0.080	15:13
03 Room Air	0.000	15:14
04 Std. Gas	0.080	15:14
05 Room Air	0.000	15:15
06 Std. Gas	0.080	15:15
07 Room Air	0.000	15:16

Lot No = 17316080A2
Cyl No = 28
Exp Date = 08/05/2018
County = 08 Oper No. = 666666

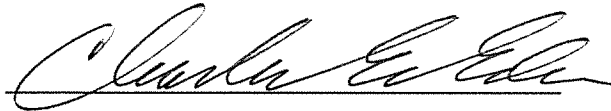

Operator Signature
CHARLES EDER

Remarks:

CALIBRATION CHECK
0.080 AC