



**ETHANOL GAS STANDARD CYLINDER REPORT**  
 OFFICE OF ATTORNEY GENERAL  
 CRIME LABORATORY DIVISION  
 Toxicology Section/Breath Alcohol Program  
 SFN 59282 (10-2013)

Chemical Test Operator Name (Print) <i>Darcy Aberle</i>		
Location <i>Williston LEC</i>		Intoxilyzer® Serial Number <i>80-005946</i>
Gas Lot Number <i>03417080A3</i>	Gas Cylinder Number <i>34</i>	Gas Expiration Date <i>03/05/2019</i>

Check When Done:

1. Scan/Enter Gas Cylinder Information (Level 1, Function S).

2. Perform an ACA Test (Level 1, Function C).

*Box was checked by RSN @ Crime Lab as ACA test was performed* Write in the ACA Test Results (Reported to 3 Digits, Ex. 0.081).

A. *0.082* AC

B. *0.081* AC

C. *0.081* AC

3. Complete the Top Portion of the Intoxilyzer® Record (SFN50496, Form 120-G) and Place it by the Intoxilyzer® for Use.

4. File Previous Intoxilyzer® Record (SFN50496, Form 120-G) at the Intoxilyzer® Location at the Agency.

5. Send the Following to the Crime Laboratory:

A. Ethanol Gas Standard Cylinder Report (SFN59282, Form 105-G).

B. ACA Test Record.

*Box was checked by RSN at Crime Lab as Form 105-G and ACA Test were received.*

Chemical Test Operator Signature <i>Darcy Aberle</i>	Date <i>10-30-17</i>
Reviewed By (Crime Laboratory Use Only) <i>Roberta Juarez-Nimmo</i>	Date <i>10-31-17</i>

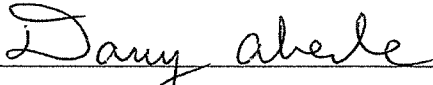
Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-005946  
Location = WILL      8164.14.00 09/16  
10/30/2017      21:47

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	21:48
02 Std. Gas	0.082	21:48
03 Room Air	0.000	21:48
04 Std. Gas	0.081	21:49
05 Room Air	0.000	21:49
06 Std. Gas	0.081	21:50
07 Room Air	0.000	21:50

Lot No = 03417080A3  
Cyl No = 34  
Exp Date = 03/05/2019  
County = 53      Oper No. = 010101

  
Operator Signature  
DARCY ABERLE

Remarks: