



INTOXILYZER® 8000 INSTALLATION AND REPAIR CHECKOUT
 OFFICE OF ATTORNEY GENERAL
 CRIME LABORATORY DIVISION
 Toxicology Section/Breath Alcohol Program
 SFN 59281 (03-2014)

Serial Number <i>80-005945</i>	Instrument Location <i>New Town New T</i>
Reason for Install/Repair <input type="checkbox"/> Install After Receiving From Crime Laboratory <input checked="" type="checkbox"/> Install After Location Change <input type="checkbox"/> Other (Specify) _____	

Check When Done:

- 1. Surge Protector Installed/Properly Grounded.
- 2. Telephone Line Connected to Intoxilyzer® 8000.
- 3. Breath Tube Heated.
- 4. Enter Preliminary Data (i.e. Date, Time, DST (Y), and Location; Level 2, Function E).
- 5. Scan/Enter Gas Cylinder Information (Level 1, Function S).
- 6. Run Tests:
 - A. Print Test (Level 1, Function P).
 - B. ACA Test (Level 1, Function C).
 - C. RFI Test (CMS Mode or Level 1, Function B or C; Key Radio During Test).
- 7. Repair and/or Maintenance Performed (if Any): *none*
- 8. Complete the Top Portion of the Intoxilyzer® Record (SFN50496, Form 120-G) and Place it by the Intoxilyzer® for Use.
- 9. File Previous Intoxilyzer® Record (SFN50496, Form 120-G) at the Intoxilyzer® Location at the Agency.
- 10. Send the Following to the Crime Laboratory: Completed Intoxilyzer®8000 Installation and Repair Checkout (SFN59281, Form 104-G), Print Test, ACA Test, and RFI Test.

*Documents were received at the Crime laboratory 2/5/18 RBN.
 Form 120-G was also sent in, but not needed for install. RBN*

Field Inspector Signature <i>[Signature]</i>	Date <i>1-31-18</i>
Reviewed By (Crime Laboratory Use Only) <i>Roberta Huegu-Tummo</i>	Date <i>2-7-18</i>

This record on file at the Office of Attorney General, Crime Laboratory Division, in the County of Burleigh, North Dakota, is certified to be a true and correct copy of the documents received.
 Certified By: *Roberta Huegu-Tummo* Certified Date: *2-7-18*

Rec'd NDCL 2/5/18 dd


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-005945
Location = NEWT 8164.14.00 09/16
01/31/2018 17:36

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	17:37
02 Std. Gas	0.081	17:37
03 Room Air	0.000	17:38
04 Std. Gas	0.081	17:38
05 Room Air	0.000	17:39
06 Std. Gas	0.080	17:39
07 Room Air	0.000	17:40

Lot No = 17316080A2
Cyl No = 35
Exp Date = 08/05/2018
County = 60 Oper No. = 132000

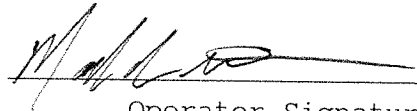
Operator Signature
MARK NEW HOLY

Remarks: *ACA Test*


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-005945
Location = NEWT 8164.14.00 09/16
01/31/2018 17:41

Test	AC	Time
01 Room Air	RFI*	17:43
02 Room Air	0.000	17:44

*Invalid Test
Inhibited - RFI

Sub Name = TES, TES T
Sub DOB = 01/01/2000
Sub Sex = Male Weight = 150
Test = DUI Cit = 123456789
Dr. Lic. = ND/TES001234
Lot No = 17316080A2
Cyl No = 35
Expiration Date = 08/05/2018
County = 60 Oper No. = 132000

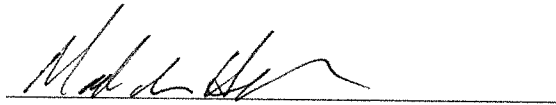

Operator Signature
MARK NEW HOLY

Remarks: *RFI Test*

