

INTOXILYZER® 8000 INSPECTION

Intoxilyzer® 8000 Serial Number: 80-005940

Location: NEWT
BIA Fort BERTHOLD
(TAT)

A. General Setup and Checks:

1. Diagnostic Tests Pass and Instrument in Ready Mode
2. Breath tube heated
3. Date, time and location code (Level 2,E). Re-set if necessary.
4. Print test (Level 1,P). Sign and attach test record.
5. Tank monitor (Level 3,D,G). Display and Regulator within 50 psi.
 - a. Display: 676 psi Regulator: 650 psi

B. Tests (Sign and attach test records):

Configure simulator for the following tests (Level 1,S).

1. Wet Calibration Check

- a. Low AC (Level 1,C): Known Value ≤ 0.03 AC: 0.020 AC
Sim. Ser #: DR3375 Lot #: 17160 Exp. Date: 6-8-19
- b. High AC (Level 1,C): Known Value ≥ 0.25 AC: 0.300 AC
Sim. Ser #: DR7111 Lot #: 17350 Exp. Date: 10-11-19

Configure dry gas standard for the following tests (Level 1,S).

2. Interferent Check (Level 1,B): Known Value: 0.10 AC + 0.05% Acetone
Sim. Ser #: DR7089 Lot #: JCS 5 Exp. Date: NA
 Display reads "Interferent Detect"
3. RFI Check (CMS Mode) Display reads "RFI Detect"
4. Dry Calibration Check (Level 1,C): Known Value 0.080 AC
Gas Cyl Lot #: 34917080A3 Cyl #: 19 Exp. Date: 2-5-20

Test 1 <u>0.079</u>	Test 4 <u>0.079</u>	Test 7 <u>0.079</u>
Test 2 <u>0.079</u>	Test 5 <u>0.079</u>	Test 8 <u>0.079</u>
Test 3 <u>0.078</u>	Test 6 <u>0.080</u>	Test 9 <u>0.079</u>
Average <u>0.079</u>		

C. Remarks/Maintenance: REPLACED SIMULATOR RETURN O-RING,
ANNUAL INSPECTION

Instrument is acceptable to be used in the field.

Breath Analyst Signature

Charles E. Esh

NA

Date

5.3.2019

NA

Reviewed by

Date

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-005940
Location = NEWT 8164.14.00 09/16
05/03/2019 12:09

***** Printer Test *****

abcdefghijklmnopqrstuvwxyz1234567890-|=|
ABCDEFGHIJKLMNPOQRSTUVWXYZ!@#\$\$%^&*()_+?

abcdefghijklmnopqrstuvwxyz1234567890-|=|
ABCDEFGHIJKLMNPOQRSTUVWXYZ!@#\$\$%^&*()_+?

Current Instrument Setup


Data Entry Mode: Enabled
Start Test Sequence: DABACABA
Display Prelim Rslt? Yes
Display Third Digit? Yes
Inhib Printer(Y/N)? No
Display Volume? No
Disable On Memfull? Yes
of Print Copies? 1
Select Std (D/W/I)? Dry
Standard Value? 0.080
Standard Lot #? 34917080A3
Standard Cyl #? 19
Standard Expiration? 02/05/2020
Oper No? 666666

Flow Cal. Date: 06/21/2013
Slope 688
Intercept -659053

IR Calibration Date: 04/02/2015

	3um	9um
0th Coef(*100):	-18917	-25088
1st Coef(*100):	267616	136158
2nd Coef(*100):	3374	1493
H2O adj(mg/l*10k):	427	417

***** Printer Test End *****

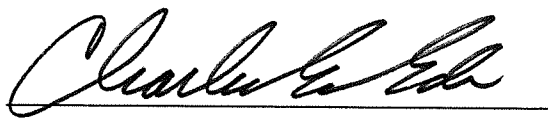

Operator Signature
CHARLES EDER

Remarks:

PRINT TEST

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-005940
Location = NEWT 8164.14.00 09/16
05/03/2019 12:11

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	12:11
02 Std. Sol.	0.020	12:12
03 Room Air	0.000	12:13
04 Std. Sol.	0.019	12:13
05 Room Air	0.000	12:14
06 Std. Sol.	0.020	12:15
07 Room Air	0.000	12:15

08 Sim Temp = 34.0°C

Simul Ser No = DR3375
Std Sol No = 17160
County = 31 Oper No. = 666666



Operator Signature
CHARLES EDER

Remarks:

Low AC
0.020 AC

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-005940
Location = NEWT 8164.14.00 09/16
05/03/2019 12:16

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	12:17
02 Std. Sol.	0.294	12:17
03 Room Air	0.000	12:18
04 Std. Sol.	0.296	12:19
05 Room Air	0.000	12:19
06 Std. Sol.	0.297	12:20
07 Room Air	0.000	12:20

08 Sim Temp = 34.0°C

Simul Ser No = DR7111

Std Sol No = 17350

County = 31

Oper No. = 666666



Operator Signature

CHARLES EDER

Remarks:

HIGH AC

0.300 AC

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-005940
Location = NEWT 8164.14.00 09/16
05/03/2019 11:59

Test	AC	Time
01 Room Air	0.000	12:00
02 *Subject Test	INT*	12:01
03 Room Air	0.000	12:01

*Invalid Test
Interferent Detected

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 01/01/1982
Sub Sex = Female Weight = 150
Test = DUI Cit = INTERFERENT
Dr. Lic. = ND/DIS821456
Lot No = 34917080A3
Cyl No = 19
Expiration Date = 02/05/2020
County = 31 Oper No. = 666666



Operator Signature
CHARLES EDER

Remarks: *Interferent Check*
0.10AC + 0.05% ACETONE

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-005940
Location = NEWT 8164.14.00 09/16
05/03/2019 12:02

Test	AC	Time
01 Diagnostic	OK	12:03
02 Room Air	0.000	12:04
03 *Subject Test	RFI*	12:04
04 Room Air	0.000	12:04

*Invalid Test
 Inhibited - RFI

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 01/01/1982
Sub Sex = Female Weight = 150
Test = DUI Cit = RFI CHECK
Dr. Lic. = ND/DIS821456
Lot No = 34917080A3
Cyl No = 19
Expiration Date = 02/05/2020
County = 31 Oper No. = 666666

I followed the Approved Method and the instructions displayed by the Intoxilyzer in conducting this test.



Operator Signature
CHARLES EDER

Remarks:

RFI CHECK

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-005940
Location = NEWT 8164.14.00 09/16
05/03/2019 12:05

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	12:06
02 Std. Gas	0.079	12:06
03 Room Air	0.000	12:07
04 Std. Gas	0.079	12:07
05 Room Air	0.000	12:07
06 Std. Gas	0.078	12:08
07 Room Air	0.000	12:08

Lot No = 34917080A3
Cyl No = 19
Exp Date = 02/05/2020
County = 31 Oper No. = 666666



Operator Signature
CHARLES EDER

Remarks: CALIBRATION CHECK
0.080 AC

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-005940
Location = NEWT 8164.14.00 09/16
05/03/2019 12:24

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	12:24
02 Std. Gas	0.079	12:25
03 Room Air	0.000	12:25
04 Std. Gas	0.079	12:26
05 Room Air	0.000	12:26
06 Std. Gas	0.080	12:26
07 Room Air	0.000	12:27

Lot No = 34917080A3
Cyl No = 19
Exp Date = 02/05/2020
County = 31

Oper No. = 666666



Operator Signature
CHARLES EDER

Remarks:

CALIBRATION CHECK
0.080 AC

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-005940
Location = NEWT 8164.14.00 09/16
05/03/2019 12:27

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	12:28
02 Std. Gas	0.079	12:28
03 Room Air	0.000	12:29
04 Std. Gas	0.079	12:29
05 Room Air	0.000	12:29
06 Std. Gas	0.079	12:30
07 Room Air	0.000	12:30

Lot No = 34917080A3
Cyl No = 19
Exp Date = 02/05/2020
County = 31 Oper No. = 666666



Operator Signature
CHARLES EDER

Remarks: CALIBRATION CHECK
0.080 AC

Form 106-I8000