

INTOXILYZER® 8000 INSPECTION

Intoxilyzer® 8000 Serial Number: 80-006497

Location: GRAND FORKS AFB.

(GFAF)

A. General Setup and Checks:

1. Diagnostic Tests Pass and Instrument in Ready Mode
2. Breath tube heated
3. Date, time and location code (Level 2,E). Re-set if necessary.
4. Print test (Level 1,P). Sign and attach test record.
5. Tank monitor (Level 3,D,G). Display and Regulator within 50 psi.
 - a. Display: 994 psi Regulator: 1000 psi

B. Tests (Sign and attach test records):

Configure simulator for the following tests (Level 1,S).

1. Wet Calibration Check
 - a. Low AC (Level 1,C): Known Value ≤ 0.03 AC: 0.020 AC
Sim. Ser #: DR5131 Lot #: 201810D Exp. Date: 10.24.20
 - b. High AC (Level 1,C): Known Value ≥ 0.25 AC: 0.300 AC
Sim. Ser #: DR7111 Lot #: 17350 Exp. Date: 10.11.19

Configure dry gas standard for the following tests (Level 1,S).

2. Interferent Check (Level 1,B): Known Value: 0.10 AC + 0.05% Acetone
Sim. Ser #: DR7089 Lot #: ICS 5 Exp. Date: NA
 Display reads "Interferent Detect"
3. RFI Check (CMS Mode) Display reads "RFI Detect"
4. Dry Calibration Check (Level 1,C): Known Value 0.080 AC
Gas Cyl Lot #: 34418080A2 Cyl #: 51 Exp. Date: 2.5.21

Test 1	<u>0.082</u>	Test 4	<u>0.082</u>	Test 7	<u>0.083</u>
Test 2	<u>0.082</u>	Test 5	<u>0.082</u>	Test 8	<u>0.083</u>
Test 3	<u>0.082</u>	Test 6	<u>0.082</u>	Test 9	<u>0.083</u>
Average	<u>0.082</u>				

C. Remarks/Maintenance: ANNUAL INSPECTION, REPLACED BREATH HOSE O-RINGS

REPLACED GAS CYLINDER, VALVE WAS BAD ON GAS CYLINDER LOT# 13518080A6, CYL: 029, EXP. 8.5.2020

Instrument is acceptable to be used in the field.

Charles E. Ehr
Breath Analyst Signature

5.15.19
Date

NA
Reviewed by

NA
Date

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-006497
Location = GFAF 8164.14.00 09/16
05/15/2019 08:48

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	08:49
02 Std. Sol.	0.019	08:50
03 Room Air	0.000	08:50
04 Std. Sol.	0.019	08:51
05 Room Air	0.000	08:52
06 Std. Sol.	0.020	08:52
07 Room Air	0.000	08:53

08 Sim Temp = 34.0°C

Simul Ser No = DR5131
Std Sol No = 201810D
County = 65 Oper No. = 666666



Operator Signature
CHARLES EDER

Remarks:

Low AC
0.020 AC

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-006497
Location = GFAF 8164.14.00 09/16
05/15/2019 08:43

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	08:44
02 Std. Sol.	0.305	08:45
03 Room Air	0.000	08:45
04 Std. Sol.	0.305	08:46
05 Room Air	0.000	08:47
06 Std. Sol.	0.306	08:47
07 Room Air	0.000	08:48

08 Sim Temp = 34.0°C

Simul Ser No = DR7111

Std Sol No = 17350

County = 65

Oper No. = 666666



Operator Signature

CHARLES EDER

Remarks:

HIGH AC
0.300 AC

Form 106-I8000


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-006497
Location = GFAF 8164.14.00 09/16
05/15/2019 08:34

Test	AC	Time
01 Room Air	0.000	08:35
02 *Subject Test	INT*	08:35
03 Room Air	0.000	08:36

*Invalid Test
Interferent Detected

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 01/01/1982
Sub Sex = Female Weight = 150
Test = DUI Cit = INTERFERENT
Dr. Lic. = ND/DIS821456
Lot No = 34418080A2
Cyl No = 51
Expiration Date = 02/05/2021
County = 65 Oper No. = 666666

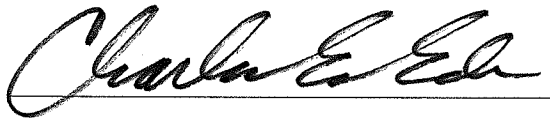

Operator Signature
CHARLES EDER

Remarks:

Interferent Check
0.10AC + 0.05% ACETONE

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-006497
Location = GFAF 8164.14.00 09/16
05/15/2019 08:37

Test	AC	Time
01 Diagnostic	OK	08:38
02 Room Air	0.000	08:38
03 *Subject Test	RFI*	08:39
04 Room Air	0.000	08:39

*Invalid Test
 Inhibited - RFI

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 01/01/1982
Sub Sex = Female Weight = 150
Test = DUI Cit = RFI CHECK
Dr. Lic. = ND/DIS821456
Lot No = 34418080A2
Cyl No = 51
Expiration Date = 02/05/2021
County = 65 Oper No. = 666666

I followed the Approved Method and the instructions displayed by the Intoxilyzer in conducting this test.



Operator Signature
CHARLES EDER

Remarks:

RFI CHECK

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-006497
Location = GFAP 8164.14.00 09/16
05/15/2019 08:30

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	08:30
02 Std. Gas	0.082	08:31
03 Room Air	0.000	08:31
04 Std. Gas	0.082	08:31
05 Room Air	0.000	08:32
06 Std. Gas	0.082	08:32
07 Room Air	0.000	08:33

Lot No = 34418080A2
Cyl No = 51
Exp Date = 02/05/2021
County = 65 Oper No. = 666666



Operator Signature
CHARLES EDER

Remarks: CALIBRATION CHECK
0.080 AC

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-006497
Location = GFAP 8164.14.00 09/16
05/15/2019 08:39

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	08:40
02 Std. Gas	0.082	08:40
03 Room Air	0.000	08:41
04 Std. Gas	0.082	08:41
05 Room Air	0.000	08:42
06 Std. Gas	0.082	08:42
07 Room Air	0.000	08:43

Lot No = 34418080A2
Cyl No = 51
Exp Date = 02/05/2021
County = 65 Oper No. = 666666



Operator Signature
CHARLES EDER

Remarks: *CALIBRATION CHECK*
0.080 AC

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-006497
Location = GFAF 8164.14.00 09/16
05/15/2019 08:54

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	08:54
02 Std. Gas	0.083	08:54
03 Room Air	0.000	08:55
04 Std. Gas	0.083	08:55
05 Room Air	0.000	08:56
06 Std. Gas	0.083	08:56
07 Room Air	0.000	08:57

Lot No = 34418080A2
Cyl No = 51
Exp Date = 02/05/2021
County = 65

Oper No. = 666666



Operator Signature
CHARLES EDER

Remarks:

CALIBRATION CHECK
0.080 AC

Form 106-I8000