



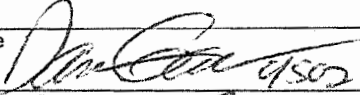

ETHANOL GAS STANDARD CYLINDER REPORT

NORTH DAKOTA OFFICE OF ATTORNEY GENERAL
CRIME LABORATORY DIVISION-TOXICOLOGY SECTION/BREATH ALCOHOL PROGRAM
SFN 59282 (08/2018)

Chemical Test Operator Name (Print) Daniel Garvens		
Location CAVA - Cavalier, ND		Intoxilyzer® Serial Number 80-005937
Gas Lot Number 34418080A2	Gas Cylinder Number 057	Gas Expiration Date 02/05/2021

Check When Done:

- 1. Scan/Enter Gas Cylinder Information (Level 1, Function S).
- 2. Perform an ACA Test (Level 1, Function C).
Write in the ACA Test Results (Reported to 3 Digits, Ex. 0.081).
 - A. 0.025 AC
 - B. 0.025 AC
 - C. 0.025 AC
- 3. Complete the Top Portion of the Intoxilyzer® Record (SFN50496, Form 120-G) and Place it by the Intoxilyzer® for Use.
- 4. File Previous Intoxilyzer® Record (SFN50496, Form 120-G) at the Intoxilyzer® Location at the Agency.
- 5. Send the Following to the Crime Laboratory:
 - A. Ethanol Gas Standard Cylinder Report (SFN59282, Form 105-G).
 - B. ACA Test Record.

Chemical Test Operator Signature 	Date 9/3/2019
Reviewed By (Crime Laboratory Use Only) 	Date 9/3/2019

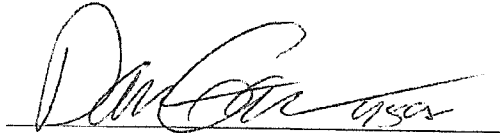
Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-005937
Location = CAVA 8164.14.00 09/16
09/03/2019 07:29

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	07:29
02 Std. Gas	0.078	07:30
03 Room Air	0.000	07:30
04 Std. Gas	0.078	07:30
05 Room Air	0.000	07:31
06 Std. Gas	0.078	07:31
07 Room Air	0.000	07:32

Lot No = 34418080A2
Cyl No = 57
Exp Date = 02/05/2021
County = 34 Oper No. = 060849



Operator Signature
DANIEL C. GARVENS

Remarks: *New Gas Cylinder Install*

Form 106-I8000