



# ETHANOL GAS STANDARD CYLINDER REPORT

NORTH DAKOTA OFFICE OF ATTORNEY GENERAL

CRIME LABORATORY DIVISION-TOXICOLOGY SECTION/BREATH ALCOHOL PROGRAM

SFN 59282 (08/2018)

Chemical Test Operator Name (Print) <b>Brian Mehlhoff</b>		
Location <b>Divide County Sheriff</b>	Intoxilyzer® Serial Number <b>80-005957</b>	
Gas Lot Number <b>2411908 ACE</b>	Gas Cylinder Number <b>32</b>	Gas Expiration Date <b>11-5-2021</b>

*OAI 1.30.20*

Check When Done:

- 1. Scan/Enter Gas Cylinder Information (Level 1, Function S).
- 2. Perform an ACA Test (Level 1, Function C). **ACA TEST (DRY CAL CHECK) SUBMITTED**  
 Write in the ACA Test Results (Reported to 3 Digits, Ex. 0.081).  
*CEE*  
*1.30.20*  
 A. 0.080 AC  
 B. 0.079 AC  
 C. 0.079 AC
- 3. Complete the Top Portion of the Intoxilyzer® Record (SFN50496, Form 120-G) and Place it by the Intoxilyzer® for Use.
- 4. File Previous Intoxilyzer® Record (SFN50496, Form 120-G) at the Intoxilyzer® Location at the Agency.
- 5. Send the Following to the Crime Laboratory:
  - A. Ethanol Gas Standard Cylinder Report (SFN59282, Form 105-G).
  - B. ACA Test Record.

Chemical Test Operator Signature <b>Brian Mehlhoff</b>	Date <b>1-27-20</b>
Reviewed By (Crime Laboratory Use Only) <b>Charles E. Ehn</b>	Date <b>1.30.20</b>


Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-005957  
Location = CROS      8164.14.00 09/16  
01/27/2020      13:27

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	13:28
02 Std. Gas	0.080	13:28
03 Room Air	0.000	13:29
04 Std. Gas	0.079	13:29
05 Room Air	0.000	13:29
06 Std. Gas	0.079	13:30
07 Room Air	0.000	13:30

Lot No = 24119080A1  
Cyl No = 32  
Exp Date = 11/05/2021  
County = 12      Oper No. = 050710

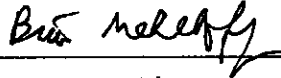

Operator Signature  
BRIAN MEHLHOFF

Remarks:

Form 106-I8000