

INTOXILYZER® 8000 INSPECTION

Intoxilyzer® 8000 Serial Number: 80-006490

Location: TOXL

A. General Setup and Checks:

1.  Diagnostic Tests Pass and Instrument in Ready Mode
2.  Breath tube heated
3.  Date, time and location code (Level 2,E). Re-set if necessary.
4.  Print test (Level 1,P). Sign and attach test record.
5.  Tank monitor (Level 3,D,G). Display and Regulator within 50 psi.
  - a. Display: 515 psi Regulator: 500 psi

B. Tests (Sign and attach test records):

Configure simulator for the following tests (Level 1,S).

1. Wet Calibration Check

- a.  Low AC (Level 1,C): Known Value  $\leq 0.03$  AC: 0.020 AC  
Sim. Ser #: MP3064 Lot #: 20070 Exp. Date: 2.13.22
- b.  High AC (Level 1,C): Known Value  $\geq 0.25$  AC: 0.300 AC 0.250 AC  
Sim. Ser #: MP3067 Lot #: 201911B Exp. Date: 11.5.21

Configure dry gas standard for the following tests (Level 1,S).

2.  Interferent Check (Level 1,B): Known Value: 0.10 AC + 0.05% Acetone  
Sim. Ser #: DR3451 Lot #: ICS 6 Exp. Date: NA  
 Display reads "Interferent Detect"
3.  RFI Check (CMS Mode)  Display reads "RFI Detect"
4.  Dry Calibration Check (Level 1,C): Known Value 0.080 AC  
Gas Cyl Lot #: 24119080A1 Cyl #: 9 Exp. Date: 11.5.21

Test 1	<u>0.081</u>	Test 4	<u>0.081</u>	Test 7	<u>0.082</u>
Test 2	<u>0.081</u>	Test 5	<u>0.082</u>	Test 8	<u>0.081</u>
Test 3	<u>0.082</u>	Test 6	<u>0.082</u>	Test 9	<u>0.081</u>
Average	<u>0.081</u>				

C. Remarks/Maintenance: ANNUAL INSPECTION

Instrument is acceptable to be used in the field.

Charles E. Ed  
Breath Analyst Signature

6.9.20  
Date

NA  
Reviewed by

NA  
Date

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer            Alcohol Analyzer  
North Dakota Model 8000        SN 80-006490  
Location = TOXL                8164.14.00 09/16  
06/09/2020                      14:30

\*\*\*\*\* Printer Test \*\*\*\*\*

abcdefghijklmnopqrstuvwxy1234567890-|=|  
ABCDEFGHIJKLMNopqrstuvwxyz!@#%&^\*()\_+?

abcdefghijklmnopqrstuvwxy1234567890-|=|  
ABCDEFGHIJKLMNopqrstuvwxyz!@#%&^\*()\_+?

Current Instrument Setup


Data Entry Mode:            Enabled  
Start Test Sequence:        DABACABA  
Display Prelim Rslt?        Yes  
Display Third Digit?        Yes  
Inhib Printer(Y/N)?        No  
Display Volume?            No  
Disable On Memfull?        Yes  
# of Print Copies?         1  
Select Std (D/W/I)?        Dry  
Standard Value?            0.080  
Standard Lot #?            24119080A1  
Standard Cyl #?            9  
Standard Expiration?       11/05/2021  
Oper No?                    666666

Flow Cal. Date:            08/18/2015  
Slope                      675  
Intercept                  -635799

IR Calibration Date:       08/18/2015

	3um	9um
0th Coef(*100):	-18073	-15729
1st Coef(*100):	265337	133887
2nd Coef(*100):	3965	1560
H2O adj(mg/l*10k):	767	457

\*\*\*\*\* Printer Test End \*\*\*\*\*

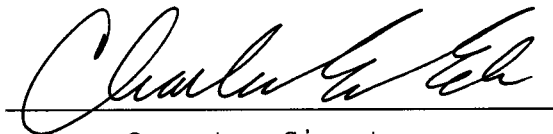

Operator Signature  
CHARLES EDER

Remarks:

PRINT TEST

Form 106-I8000

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer            Alcohol Analyzer  
North Dakota Model 8000        SN 80-006490  
Location = TOXL                    8164.14.00 09/16  
06/09/2020                            14:47

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	14:48
02 Std. Sol.	0.019	14:49
03 Room Air	0.000	14:49
04 Std. Sol.	0.018	14:50
05 Room Air	0.000	14:51
06 Std. Sol.	0.019	14:51
07 Room Air	0.000	14:52

08 Sim Temp = 34.0°C

Simul Ser No = MP3064  
Std Sol No = 20070  
County = 08                            Oper No. = 666666



Operator Signature  
CHARLES EDER

Remarks:    *LOW AC*  
                  *0.020 AC*

Form 106-I8000

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-006490  
Location = TOXL      8164.14.00 09/16  
06/09/2020      15:03

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	15:03
02 Std. Sol.	0.255	15:04
03 Room Air	0.000	15:05
04 Std. Sol.	0.254	15:05
05 Room Air	0.000	15:06
06 Std. Sol.	0.254	15:07
07 Room Air	0.000	15:07

08 Sim Temp = 34.0°C

Simul Ser No = MP3067  
Std Sol No = 201911B  
County = 08      Oper No. = 666666



Operator Signature  
CHARLES EDER

Remarks:      HIGH AC  
                    0.250 AC

Form 106-I8000


Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-006490  
Location = TOXL      8164.14.00 09/16  
06/09/2020      15:09

Test	AC	Time
01 Room Air	0.000	15:10
02 *Subject Test	INT*	15:10
03 Room Air	0.000	15:11

\*Invalid Test  
Interferent Detected

Sub Name = DISCOVER, THE SPIRIT  
Sub DOB = 02/01/1992  
Sub Sex = Female      Weight = 150  
Test = DUI      Cit = INTERFERENT  
Dr. Lic. = ND/DIS921456  
Lot No = 24119080A1  
Cyl No = 9  
Expiration Date = 11/05/2021  
County = 08      Oper No. = 666666

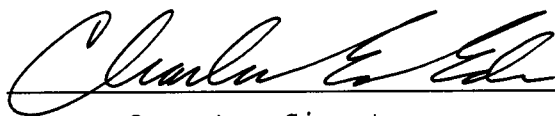

Operator Signature  
CHARLES EDER

Remarks:

*Interferent Check*  
*0.10AC + 0.05% ACETONE*

Form 106-I8000

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501


CMI, Inc. Intoxilyzer            Alcohol Analyzer  
North Dakota Model 8000        SN 80-006490  
Location = TOXL                8164.14.00 09/16  
06/09/2020                      15:11

Test	AC	Time
01 Diagnostic	OK	15:12
02 Room Air	0.000	15:13
03 *Subject Test	RFI*	15:13
04 Room Air	0.000	15:13

\*Invalid Test  
Inhibited - RFI

Sub Name = DISCOVER, THE SPIRIT  
Sub DOB = 02/01/1992  
Sub Sex = Female                Weight = 150  
Test = DUI                        Cit = RFI CHECK  
Dr. Lic. = ND/DIS921456  
Lot No = 24119080A1  
Cyl No = 9  
Expiration Date = 11/05/2021  
County = 08                      Oper No. = 666666

I followed the Approved Method and the instructions displayed by the Intoxilyzer in conducting this test.

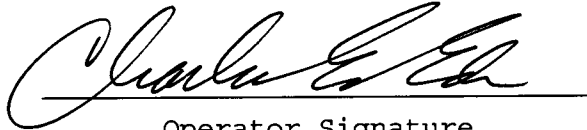

Operator Signature  
CHARLES EDER

Remarks:

*RFI CHECK*

Form 106-I8000


Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-006490  
Location = TOXL      8164.14.00 09/16  
06/09/2020      15:14

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	15:14
02 Std. Gas	0.081	15:14
03 Room Air	0.000	15:15
04 Std. Gas	0.081	15:15
05 Room Air	0.000	15:16
06 Std. Gas	0.082	15:16
07 Room Air	0.000	15:17

Lot No = 24119080A1  
Cyl No = 9  
Exp Date = 11/05/2021  
County = 08      Oper No. = 666666

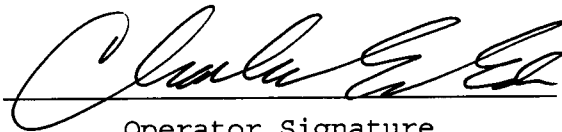

Operator Signature  
CHARLES EDER

Remarks: CALIBRATION CHECK  
0-080 AC

Form 106-I8000


Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-006490  
Location = TOXL      8164.14.00 09/16  
06/09/2020      15:17

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	15:18
02 Std. Gas	0.081	15:18
03 Room Air	0.000	15:19
04 Std. Gas	0.082	15:19
05 Room Air	0.000	15:20
06 Std. Gas	0.082	15:20
07 Room Air	0.000	15:20

Lot No = 24119080A1  
Cyl No = 9  
Exp Date = 11/05/2021  
County = 08      Oper No. = 666666

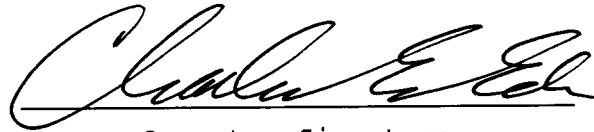

Operator Signature  
CHARLES EDER

Remarks: CALIBRATION CHECK  
0.080 AC

Form 106-I8000



Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-006490  
Location = TOXL      8164.14.00 09/16  
06/09/2020      15:21

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	15:21
02 Std. Gas	0.082	15:22
03 Room Air	0.000	15:22
04 Std. Gas	0.081	15:23
05 Room Air	0.000	15:23
06 Std. Gas	0.081	15:24
07 Room Air	0.000	15:24

Lot No = 24119080A1  
Cyl No = 9  
Exp Date = 11/05/2021  
County = 08      Oper No. = 666666



Operator Signature  
CHARLES EDER

Remarks: CALIBRATION CHECK  
0.080 AC

Form 106-I8000