

INTOXILYZER® 8000 INSPECTION

Intoxilyzer® 8000 Serial Number: 80-005948 Location: TOXL

A. General Setup and Checks:

- 1.  Diagnostic Tests Pass and Instrument in Ready Mode
- 2.  Breath tube heated
- 3.  Date, time and location code (Level 2,E). Re-set if necessary.
- 4.  Print test (Level 1,P). Sign and attach test record.
- 5.  Tank monitor (Level 3,D,G). Display and Regulator within 50 psi.
  - a. Display: 936 psi Regulator: 950psi

B. Tests (Sign and attach test records):

Configure simulator for the following tests (Level 1,S).

- 1. Wet Calibration Check
  - a.  Low AC (Level 1,C): Known Value  $\leq 0.03$  AC: 0.020 AC  
Sim. Ser #: MP3064 Lot #: 20070 Exp. Date: 2.13.22
  - b.  High AC (Level 1,C): Known Value  $\geq 0.25$  AC: 0.300 AC  
Sim. Ser #: MP5317 Lot #: 202012A Exp. Date: 12.8.22

Configure dry gas standard for the following tests (Level 1,S).

- 2.  Interferent Check (Level 1,B): Known Value: 0.10 AC + 0.05% Acetone  
Sim. Ser #: DR5114 Lot #: IC56 Exp. Date: NA  
 Display reads "Interferent Detect"
- 3.  RFI Check (CMS Mode)  Display reads "RFI Detect"
- 4.  Dry Calibration Check (Level 1,C): Known Value 0.080 AC  
Gas Cyl Lot #: 0262108041 Cyl #: 26 Exp. Date: 3.5.23

Test 1	<u>0.081</u>	Test 4	<u>0.080</u>	Test 7	<u>0.081</u>
Test 2	<u>0.080</u>	Test 5	<u>0.080</u>	Test 8	<u>0.080</u>
Test 3	<u>0.079</u>	Test 6	<u>0.080</u>	Test 9	<u>0.080</u>
Average	<u>0.080</u>				

C. Remarks/Maintenance: ANNUAL INSPECTION

Instrument is acceptable to be used in the field.

Charles E. Er  
Breath Analyst Signature

4.22.21  
Date

Katlynn Hib  
Reviewed by

4.27.21  
Date

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer            Alcohol Analyzer  
North Dakota Model 8000            SN 80-005948  
Location = TOXL                    8164.14.00 09/16  
04/22/2021                            13:11

\*\*\*\*\* Printer Test \*\*\*\*\*

abcdefghijklmnopqrstuvwxyz1234567890- = |  
ABCDEFGHIJKLMNopqrstuvwxyz!@#%^&\*()\_+?

abcdefghijklmnopqrstuvwxyz1234567890- = |  
ABCDEFGHIJKLMNopqrstuvwxyz!@#%^&\*()\_+?

Current Instrument Setup

Data Entry Mode:                    Enabled  
Start Test Sequence:                DABACABA  
Display Prelim Rslt?                Yes  
Display Third Digit?                Yes  
Inhib Printer(Y/N)?                No  
Display Volume?                    No  
Disable On Memfull?                Yes  
# of Print Copies?                  1  
Select Std (D/W/I)?                Dry  
Standard Value?                    0.080  
Standard Lot #?                    02621080A1  
Standard Cyl #?                    26  
Standard Expiration?               03/05/2023  
Oper No?                            666666

Flow Cal. Date:                    09/12/2019  
Slope                                462  
Intercept                            435439

IR Calibration Date:               09/12/2019  
    3um                            9um

-----  
0th Coef(\*100):                    -18342                    -17009  
1st Coef(\*100):                    275562                    137785  
2nd Coef(\*100):                    2253                        1137  
H2O adj(mg/l\*10k):                428                         318

\*\*\*\*\* Printer Test End \*\*\*\*\*



Operator Signature  
CHARLES EDER

Remarks:

PRINT TEST

Form 106-I8000

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501


CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-005948  
Location = TOXL      8164.14.00 09/16  
04/22/2021      13:14

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	13:15
02 Std. Sol.	0.018	13:15
03 Room Air	0.000	13:16
04 Std. Sol.	0.019	13:17
05 Room Air	0.000	13:17
06 Std. Sol.	0.020	13:18
07 Room Air	0.000	13:18

08 Sim Temp = 34.0°C

Simul Ser No = MP3064  
Std Sol No = 20070  
County = 08      Oper No. = 666666

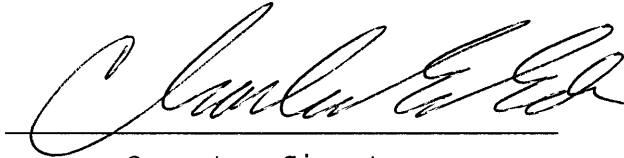

Operator Signature  
CHARLES EDER

Remarks:

Low AC  
0.020 AC

Form 106-I8000

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-005948  
Location = TOXL      8164.14.00 09/16  
04/22/2021      13:19

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	13:20
02 Std. Sol.	0.298	13:20
03 Room Air	0.000	13:21
04 Std. Sol.	0.298	13:22
05 Room Air	0.000	13:22
06 Std. Sol.	0.299	13:23
07 Room Air	0.000	13:24

08 Sim Temp = 34.0°C

Simul Ser No = MP5317  
Std Sol No = 202012A  
County = 08      Oper No. = 666666



Operator Signature  
CHARLES EDER

Remarks:

HIGH AC  
0.300 AC

Form 106-I8000

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer            Alcohol Analyzer  
North Dakota Model 8000            SN 80-005948  
Location = TOXL                    8164.14.00 09/16  
04/22/2021                            13:25

Test	AC	Time
01 Room Air	0.000	13:26
02 *Subject Test	INT*	13:26
03 Room Air	0.000	13:27

\*Invalid Test  
Interferent Detected

Sub Name = DISCOVER, THE SPIRIT  
Sub DOB = 02/01/1992  
Sub Sex = Female                    Weight = 150  
Test = DUI                            Cit = INTERFERENT  
Dr. Lic. = ND/DIS921456  
Lot No = 02621080A1  
Cyl No = 26  
Expiration Date = 03/05/2023  
County = 08                            Oper No. = 666666



Operator Signature  
CHARLES EDER

Remarks:

*Interferent Check*  
*0.10AC + 0.05% ACETONE*

Form 106-I8000

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer            Alcohol Analyzer  
North Dakota Model 8000        SN 80-005948  
Location = TOXL                 8164.14.00 09/16  
04/22/2021                        13:27

Test	AC	Time
01 Diagnostic	OK	13:28
02 Room Air	0.000	13:28
03 *Subject Test	RFI*	13:28
04 Room Air	0.000	13:29

\*Invalid Test  
Inhibited - RFI

Sub Name = DISCOVER, THE SPIRIT  
Sub DOB = 02/01/1992  
Sub Sex = Female                 Weight = 150  
Test = DUI                         Cit = RFI CHECK  
Dr. Lic. = ND/DIS921456  
Lot No = 02621080A1  
Cyl No = 26  
Expiration Date = 03/05/2023  
County = 08                        Oper No. = 666666

I followed the Approved Method and the instructions displayed by the Intoxilyzer in conducting this test.



Operator Signature  
CHARLES EDER

Remarks: *RFI CHECK*

Form 106-I8000

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-005948  
Location = TOXL      8164.14.00 09/16  
04/22/2021      13:29

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	13:30
02 Std. Gas	0.081	13:30
03 Room Air	0.000	13:31
04 Std. Gas	0.080	13:31
05 Room Air	0.000	13:32
06 Std. Gas	0.079	13:32
07 Room Air	0.000	13:33

Lot No = 02621080A1  
Cyl No = 26  
Exp Date = 03/05/2023  
County = 08      Oper No. = 666666



Operator Signature  
CHARLES EDER

Remarks: CALIBRATION CHECK  
0.080 AC

Form 106-I8000


Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-005948  
Location = TOXL      8164.14.00 09/16  
04/22/2021      13:33

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	13:33
02 Std. Gas	0.080	13:34
03 Room Air	0.000	13:34
04 Std. Gas	0.080	13:34
05 Room Air	0.000	13:35
06 Std. Gas	0.080	13:35
07 Room Air	0.000	13:36

Lot No = 02621080A1  
Cyl No = 26  
Exp Date = 03/05/2023  
County = 08      Oper No. = 666666

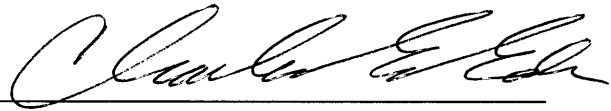

Operator Signature  
CHARLES EDER

Remarks: CALIBRATION CHECK  
0.080 AC

Form 106-I8000




Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-005948  
Location = TOXL      8164.14.00 09/16  
04/22/2021      13:36

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	13:37
02 Std. Gas	0.081	13:37
03 Room Air	0.000	13:38
04 Std. Gas	0.080	13:38
05 Room Air	0.000	13:39
06 Std. Gas	0.080	13:39
07 Room Air	0.000	13:39

Lot No = 02621080A1  
Cyl No = 26  
Exp Date = 03/05/2023  
County = 08      Oper No. = 666666

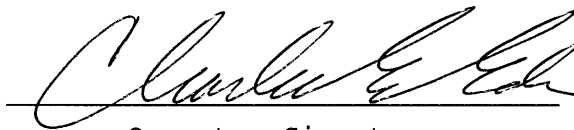

Operator Signature  
CHARLES EDER

Remarks:

CALIBRATION CHECK  
0.080 AC

Form 106-I8000