



INTOXILYZER® 8000 INSTALLATION AND REPAIR CHECKOUT
 NORTH DAKOTA OFFICE OF ATTORNEY GENERAL
 CRIME LABORATORY DIVISION-TOXICOLOGY SECTION/BREATH ALCOHOL PROGRAM
 SFN 59281 (06/2018)

Serial Number <i>SN 80-00495A</i>	Instrument Location <i>WCPJ</i>
Reason for Install/Repair <input checked="" type="checkbox"/> Install After Receiving From Crime Laboratory <input type="checkbox"/> Install After Location Change <input type="checkbox"/> Other (Specify) _____	

Check When Done:

- 1. Surge Protector Installed/Property Grounded.
- 2. Telephone Line Connected to Intoxilyzer® 8000.
- 3. Breath Tube Heated.
- 4. Enter Preliminary Data (i.e. Date, Time, DST (Y), and Location; Level 2, Function E).
- 5. Scan/Enter Gas Cylinder Information (Level 1, Function S).
- 6. Run Tests:
 - A. Print Test (Level1, Function P).
 - B. ACA Test (Level 1, Function C).
 - C. Radio Frequency Interference (RFI) Test (CMS Mode or Level 1, Function B or C; Key Radio During Test).
- 7. Repair and/or Maintenance Performed (if any): _____
- 8. Complete the Top Portion of the Intoxilyzer® Record (SFN50496, Form 120-G) and Place it by the Intoxilyzer® for Use.
- 9. File Previous Intoxilyzer® Record (SFN504096, Form 120-G) at the Intoxilyzer® Location at the Agency.
- 10. Send the Following to the Crime Laboratory: Completed Intoxilyzer® 8000 Installation and Repair Checkout (SFN59281, Form 104-G), Print Test, ACA Test, and RFI Test.

Field Inspector Signature <i>M. Johnson</i>	Date <i>05/20/22</i>
--	-------------------------

Crime Laboratory Use Only

This installation has been reviewed and the instrument is approved to be used for the analysis of breath to determine alcohol concentration from the date the Field Inspector performed the installation. This record on file at the Office of Attorney General, Crime Laboratory Division, in the County of Burleigh, North Dakota, is certified to be a true and correct copy of the documents received.

Reviewed/Certified By <i>Janelle Putschiller</i>	Certified Date <i>06 June 22</i>
---	-------------------------------------

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004954
Location = WCPD 8164.14.00 09/16
05/20/2022 11:31

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	11:32
02 Std. Gas	0.082	11:32
03 Room Air	0.000	11:33
04 Std. Gas	0.082	11:33
05 Room Air	0.000	11:34
06 Std. Gas	0.082	11:34
07 Room Air	0.000	11:34


Lot No = 33321080A4

Cyl No = 11

Exp Date = 01/05/2024

County = 27

Oper No. = 133050

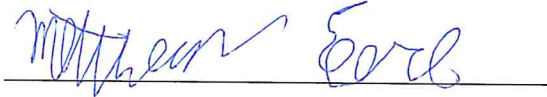

Operator Signature
MATTHEW EARL

Remarks:

IN STPV

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004954
Location = WCPD 8164.14.00 09/16
05/20/2022 11:36

Test	AC	Time
01 Room Air	0.000	11:39
02 *Subject Test	RFI*	11:39
03 Room Air	0.000	11:40

*Invalid Test
Inhibited - RFI

Sub Name = EARL, MATTHEW EMMETT
Sub DOB = 05/13/2000
Sub Sex = Male Weight = 125
Test = OTH Cit = NA
Dr. Lic. = ND/EAR004876
Lot No = 33321080A4
Cyl No = 11
Expiration Date = 01/05/2024
County = 27 Oper No. = 133050



Operator Signature
MATTHEW EARL

Remarks:

INSTALL