

**INTOXILYZER® 8000 INSPECTION**

Intoxilyzer® 8000 Serial Number: 80-004188 Location: TOXL

A. General Setup and Checks:

- 1.  Diagnostic Tests Pass and Instrument in Ready Mode
- 2.  Breath tube heated
- 3.  Date, time and location code (Level 2,E). Re-set if necessary.
- 4.  Print test (Level 1,P). Sign and attach test record.
- 5.  Tank monitor (Level 3,D,G). Display and Regulator within 50 psi.
  - a. Display: 226 psi Regulator: 280 psi

B. Tests (Sign and attach test records):

- Configure simulator for the following tests (Level 1,S).
  - 1. Wet Calibration Check
    - a.  Low AC (Level 1,C): Known Value  $\leq 0.03$  AC: 0.020 AC  
Sim. Ser #: MP3070 Lot #: 22CP41184 Exp. Date: 31/1/24
    - b.  High AC (Level 1,C): Known Value  $\geq 0.25$  AC: 0.300 AC  
Sim. Ser #: MP3068 Lot #: 2022 01 F Exp. Date: 1/18/24

Configure dry gas standard for the following tests (Level 1,S).

- 2.  Interferent Check (Level 1,B): Known Value: 0.10 AC + 0.05% Acetone  
Sim. Ser #: DR3847 Lot #: 1CS 6 Exp. Date: N/A  
 Display reads "Interferent Detect"
- 3.  RFI Check (CMS Mode)  Display reads "RFI Detect" 0.080
- 4.  Dry Calibration Check (Level 1,C): Known Value 0.0800 AC  
Gas Cyl Lot #: 02621080A1 Cyl #: 16 Exp. Date: 3/5/23

Test 1	<u>0.078</u>	Test 4	<u>0.078</u>	Test 7	<u>0.077</u>
Test 2	<u>0.077</u>	Test 5	<u>0.077</u>	Test 8	<u>0.078</u>
Test 3	<u>0.077</u>	Test 6	<u>0.077</u>	Test 9	<u>0.077</u>
Average <u>0.077</u>					

C. Remarks/Maintenance: Annual inspection.  
Replaced breath tube and simulator return o-rings.

Instrument is acceptable to be used in the field.

Janelle Pertschiller  
Breath Analyst Signature

6/22/22  
Date

[Signature]  
Reviewed by

062322  
Date



Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-004188  
Location = TOXL      8164.14.00 09/16  
06/22/2022      11:41

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	11:41
02 Std. Sol.	0.019	11:42
03 Room Air	0.000	11:43
04 Std. Sol.	0.019	11:43
05 Room Air	0.000	11:44
06 Std. Sol.	0.019	11:45
07 Room Air	0.000	11:45

08 Sim Temp = 34.0°C

Simul Ser No = MP3070  
Std Sol No = 22CP41184  
County = 08      Oper No. = 132128

  
Operator Signature  
JANELLE PORTSCHELLER

Remarks: *Low Ac Check 0.020*

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-004188  
Location = TOXL      8164.14.00 09/16  
06/22/2022      12:00

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	12:01
02 Std. Sol.	0.293	12:02
03 Room Air	0.000	12:02
04 Std. Sol.	0.296	12:03
05 Room Air	0.000	12:04
06 Std. Sol.	0.295	12:04
07 Room Air	0.000	12:05

08 Sim Temp = 34.0°C

Simul Ser No = MP3068

Std Sol No = 202201F

County = 08

Oper No. = 132128



Operator Signature  
JANELLE PORTSCHELLER

Remarks: *High Ac Check 0.300*

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer            Alcohol Analyzer  
North Dakota Model 8000        SN 80-004188  
Location = TOXL                    8164.14.00 09/16  
06/22/2022                            12:11

Test	AC	Time
01 Room Air	0.000	12:12
02 *Subject Test	INT*	12:12
03 Room Air	0.000	12:13

\*Invalid Test  
Interferent Detected

Sub Name = DISCOVER, THE SPIRIT  
Sub DOB = 01/01/1982  
Sub Sex = Female                    Weight = 150  
Test = DUI                            Cit = INTERFERENT CHK  
Dr. Lic. = ND/DIS821456  
Lot No = 02621080A1  
Cyl No = 16  
Expiration Date = 03/05/2023  
County = 08                            Oper No. = 132128



Operator Signature  
JANELLE PORTSCHELLER

Remarks: *RTMP Interferent Check*

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer            Alcohol Analyzer  
North Dakota Model 8000            SN 80-004188  
Location = TOXL                    8164.14.00 09/16  
06/22/2022                            12:13


Test	AC	Time
01 Diagnostic	OK	12:14
02 Room Air	0.000	12:15
03 *Subject Test	0.000*	12:18
04 Room Air	RFI**	12:20
05 Room Air	0.000	12:20

\*Deficient Sample - Value Printed was  
Highest Obtained

\*\*Invalid Test  
Inhibited - RFI

Sub Name = DISCOVER, THE SPIRIT  
Sub DOB = 01/01/1982  
Sub Sex = Female                    Weight = 150  
Test = DUI                            Cit = RFI CHK  
Dr. Lic. = ND/DIS821456  
Lot No = 02621080A1  
Cyl No = 16  
Expiration Date = 03/05/2023  
County = 08                            Oper No. = 132128

I followed the Approved Method and the  
instructions displayed by the Intoxilyzer  
in conducting this test.

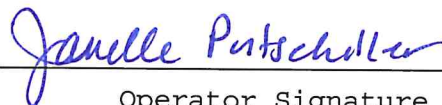

Operator Signature  
JANELLE PORTSCHELLER

Remarks: *RFI check*










Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer            Alcohol Analyzer  
North Dakota Model 8000        SN 80-004188  
Location = TOXL                    8164.14.00 09/16  
06/22/2022                            12:33

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	12:34
02 Std. Gas	0.077	12:34
03 Room Air	0.000	12:35
04 Std. Gas	0.078	12:35
05 Room Air	0.000	12:35
06 Std. Gas	0.077	12:36
07 Room Air	0.000	12:36

Lot No = 02621080A1  
Cyl No = 16  
Exp Date = 03/05/2023  
County = 08                            Oper No. = 132128

  
\_\_\_\_\_  
Operator Signature  
JANELLE PORTSCHELLER

Remarks: *Dry Calibration Check 7-9*