



ETHANOL GAS STANDARD CYLINDER REPORT

NORTH DAKOTA OFFICE OF ATTORNEY GENERAL
CRIME LABORATORY DIVISION-TOXICOLOGY SECTION/BREATH ALCOHOL PROGRAM
SFN 59282 (08/2018)

Chemical Test Operator Name (Print) <i>Taylor Wray</i>		
Location <i>Burleigh Morton Detention Center</i>	Intoxilyzer® Serial Number <i>80-004945</i>	
Gas Lot Number <i>01923080A3</i>	Gas Cylinder Number <i>055</i>	Gas Expiration Date <i>02/05/2025</i>

Check When Done:

- 1. Scan/Enter Gas Cylinder Information (Level 1, Function S).
- 2. Perform an ACA Test (Level 1, Function C).
Write in the ACA Test Results (Reported to 3 Digits, Ex. 0.081).
 - A. *0.083* AC
 - B. *0.083* AC
 - C. *0.084* AC
- 3. Complete the Top Portion of the Intoxilyzer® Record (SFN50496, Form 120-G) and Place it by the Intoxilyzer® for Use.
- 4. File Previous Intoxilyzer® Record (SFN50496, Form 120-G) at the Intoxilyzer® Location at the Agency.
- 5. Send the Following to the Crime Laboratory:
 - A. Ethanol Gas Standard Cylinder Report (SFN59282, Form 105-G).
 - B. ACA Test Record.

Chemical Test Operator Signature <i>Taylor Wray</i>	Date <i>03/02/2023</i>
Reviewed By (Crime Laboratory Use Only) <i>[Signature]</i>	Date <i>03 March 2023</i>

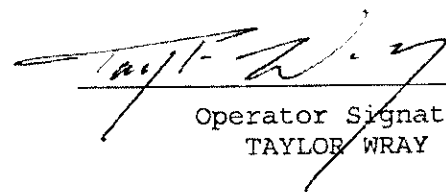
Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004945
Location = BMDC 8164.14.00 09/16
03/02/2023 09:29

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	09:30
02 Std. Gas	0.083	09:30
03 Room Air	0.000	09:31
04 Std. Gas	0.083	09:31
05 Room Air	0.000	09:31
06 Std. Gas	0.084	09:32
07 Room Air	0.000	09:32

Lot No = 01923080A3
Cyl No = 55
Exp Date = 02/05/2025
County = 08 Oper No. = 131856


Operator Signature
TAYLOR WRAY

Remarks:

Form 106-I8000