



**ETHANOL GAS STANDARD CYLINDER REPORT**  
 NORTH DAKOTA OFFICE OF ATTORNEY GENERAL  
 CRIME LABORATORY DIVISION-TOXICOLOGY SECTION/BREATH ALCOHOL PROGRAM  
 SFN 59282 (08/2018)

Chemical Test Operator Name (Print) <i>Michael P Beaver</i>		
Location <i>Wells Co. S.O.</i>	Intoxilyzer® Serial Number <i>80-006504</i>	
Gas Lot Number <sup>MS</sup> <del><i>80-006504</i></del> <i>26021080A1</i>	Gas Cylinder Number <i>30</i>	Gas Expiration Date <i>10/05/2023</i>

Check When Done:

- 1. Scan/Enter Gas Cylinder Information (Level 1, Function S).
- 2. Perform an ACA Test (Level 1, Function C).  
 Write in the ACA Test Results (Reported to 3 Digits, Ex. 0.081).  
 A. *0080* AC  
 B. *0029* AC  
 C. *0080* AC
- 3. Complete the Top Portion of the Intoxilyzer® Record (SFN50496, Form 120-G) and Place it by the Intoxilyzer® for Use.
- 4. File Previous Intoxilyzer® Record (SFN50496, Form 120-G) at the Intoxilyzer® Location at the Agency.
- 5. Send the Following to the Crime Laboratory:  
 A. Ethanol Gas Standard Cylinder Report (SFN59282, Form 105-G).  
 B. ACA Test Record.

*JP*  
*12-Apr-22*  
*See attached*

Chemical Test Operator Signature 	Date <i>02-Apr-2022</i>
Reviewed By (Crime Laboratory Use Only) <i>Janelli Putschillo</i>	Date <i>12 Apr 22</i>

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501


CMI, Inc. Intoxilyzer Alcohol Analyzer  
North Dakota Model 8000 SN 80-006504  
Location = WESO 8164.14.00 09/16  
04/02/2022 10:14

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	10:14
02 Std. Gas	0.080	10:15
03 Room Air	0.000	10:15
04 Std. Gas	0.079	10:16
05 Room Air	0.000	10:16
06 Std. Gas	0.080	10:17
07 Room Air	0.000	10:17

Lot No = 26021080A1  
Cyl No = 30  
Exp Date = 10/05/2023  
County = 52

Oper No. = 131495

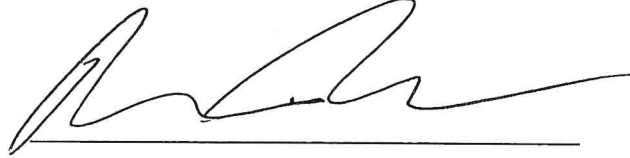

Operator Signature  
MICHAEL BEAVER

Remarks:

Form 106-I8000