



ETHANOL BREATH STANDARD CYLINDER REPORT

OFFICE OF ATTORNEY GENERAL
CRIME LABORATORY DIVISION
Toxicology Section/Breath Alcohol Program
SFN 59282 (5-2011)

Chemical Test Operator Name (Print) <i>Trevina, Joshua T</i>		
Location <i>Watford City Police Department</i>		Intoxilyzer® Serial Number <i>80-005958</i>
Gas Lot Number <i>14214080 AB</i>	Gas Cylinder Number <i>050</i>	Gas Expiration Date <i>07-05-2016</i>

Check When Done:

- 1. Scan/Enter Gas Cylinder Information (Level 1, Function S).
- 2. Perform an ACA Test (Level 1, Function C).
Write in the ACA Test Results (Reported to 3 Digits, Ex. 0.081).
 - A. 0. 080 AC
 - B. 0. 081 AC
 - C. 0. 080 AC
- 3. Complete the Top Portion of the Intoxilyzer® Record (SFN50496, Form 120-G) and Place it by the Intoxilyzer® for Use.
- 4. File Previous Intoxilyzer® Record (SFN50496, Form 120-G) at the Intoxilyzer® Location at the Agency.
- 5. Send the Following to the Crime Laboratory:
 - A. Ethanol Breath Standard Cylinder Report (SFN59282, Form 105-G).
 - B. ACA Test Record.

Chemical Test Operator Signature <i>[Signature]</i>	Date <i>12-31-2015</i>
Reviewed By (Crime Laboratory Use Only) <i>[Signature]</i>	Date <i>31 Dec 2015</i>

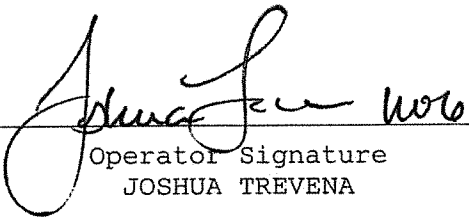
Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-005958
Location = WCPD 8164.13.00 06/09
12/31/2015 15:25

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	15:26
02 Std. Gas	0.080	15:26
03 Room Air	0.000	15:27
04 Std. Gas	0.081	15:27
05 Room Air	0.000	15:28
06 Std. Gas	0.080	15:28
07 Room Air	0.000	15:29

Lot No = 14214080A6
Cyl No = 050
Exp Date = 07/05/2016
County = 27 Oper No. = 130612


Operator Signature
JOSHUA TREVENA

Remarks: *new gas standard cylinder*