


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004204
Location = TOXL 8164.14.00 09/16
06/28/2023 12:32

Test	AC	Time
01 Room Air	0.000	12:32
02 Std. Gas	0.078	12:33
03 Room Air	0.000	12:33
04 Std. Gas	0.079	12:34
05 Room Air	0.000	12:34
06 Std. Gas	0.079	12:34
07 Room Air	0.000	12:35
08 Std. Gas	0.080	12:35
09 Room Air	0.000	12:36
10 Std. Gas	0.079	12:36
11 Room Air	0.000	12:37
12 Std. Gas	0.079	12:37
13 Room Air	0.000	12:38
14 Std. Gas	0.080	12:38
15 Room Air	0.000	12:38
16 Std. Gas	0.079	12:39
17 Room Air	0.000	12:39
18 Std. Gas	0.080	12:40
19 Room Air	0.000	12:40
20 Std. Gas	0.080	12:40
21 Room Air	0.000	12:41
22 Std. Gas	0.080	12:41
23 Room Air	0.000	12:42
24 Std. Gas	0.080	12:42
25 Room Air	0.000	12:43
26 Std. Gas	0.080	12:43
27 Room Air	0.000	12:44
28 Std. Gas	0.079	12:44
29 Room Air	0.000	12:44
30 Std. Gas	0.080	12:45
31 Room Air	0.000	12:45
32 Cal Check Stats:		
33 Average	0.0795	
34 Std Dev	0.0006	
35 Rel Std Dev	0.81%	

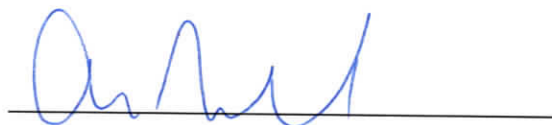

Operator Signature
ANNA NAREHOOD

Remarks: *Stability Test*