



ETHANOL GAS STANDARD CYLINDER REPORT
 NORTH DAKOTA OFFICE OF ATTORNEY GENERAL
 CRIME LABORATORY DIVISION-TOXICOLOGY SECTION/BREATH ALCOHOL PROGRAM
 SFN 59282 (08/2018)

Chemical Test Operator Name (Print) Joshua Fastnacht 368		
Location JAME Adult Upstairs	Intoxilyzer® Serial Number 80-006499	
Gas Lot Number 14323080A1	Gas Cylinder Number 24	Gas Expiration Date 6-5-2025

Check When Done:

- 1. Scan/Enter Gas Cylinder Information (Level 1, Function S).
- 2. Perform an ACA Test (Level 1, Function C).
 Write in the ACA Test Results (Reported to 3 Digits, Ex. 0.081).
 A. 0.080 AC
 B. 0.081 AC
 C. 0.081 AC
- 3. Complete the Top Portion of the Intoxilyzer® Record (SFN50496, Form 120-G) and Place it by the Intoxilyzer® for Use.
- 4. File Previous Intoxilyzer® Record (SFN50496, Form 120-G) at the Intoxilyzer® Location at the Agency.
- 5. Send the Following to the Crime Laboratory:
 - A. Ethanol Gas Standard Cylinder Report (SFN59282, Form 105-G).
 - B. ACA Test Record.

Chemical Test Operator Signature 	Date 9/21/2023
Reviewed By (Crime Laboratory Use Only) 	Date 22Sep2023

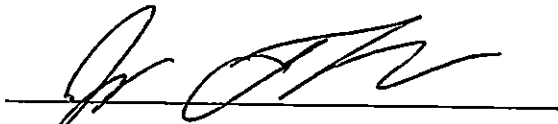
Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-006499
Location = JAME 8164.14.00 09/16
09/21/2023 04:19

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	04:19
02 Std. Gas	0.080	04:19
03 Room Air	0.000	04:20
04 Std. Gas	0.081	04:20
05 Room Air	0.000	04:21
06 Std. Gas	0.081	04:21
07 Room Air	0.000	04:22

Lot No = 14323080A1
Cyl No = 24
Exp Date = 06/05/2025
County = 47 Oper No. = 132539


Operator Signature
JOSHUA FASTNACHT

Remarks:

Tank Install

Form 106-I8000