



ETHANOL GAS STANDARD CYLINDER REPORT

NORTH DAKOTA OFFICE OF ATTORNEY GENERAL

CRIME LABORATORY DIVISION-TOXICOLOGY SECTION/BREATH ALCOHOL PROGRAM

SFN 59282 (08/2018)

Chemical Test Operator Name (Print) <i>Cory Mortensen</i>		
Location <i>Cavalier</i>	Intoxilyzer® Serial Number <i>80-006494</i>	
Gas Lot Number <i>14323080A4</i>	Gas Cylinder Number <i>6</i>	Gas Expiration Date <i>06/05/2025</i>

Check When Done:

- 1. Scan/Enter Gas Cylinder Information (Level 1, Function S).
- 2. Perform an ACA Test (Level 1, Function C).
Write in the ACA Test Results (Reported to 3 Digits, Ex. 0.081).
A. 0. *083* AC
B. 0. *083* AC
C. 0. *083* AC
- 3. Complete the Top Portion of the Intoxilyzer® Record (SFN50496, Form 120-G) and Place it by the Intoxilyzer® for Use.
- 4. File Previous Intoxilyzer® Record (SFN50496, Form 120-G) at the Intoxilyzer® Location at the Agency.
- 5. Send the Following to the Crime Laboratory:
 - A. Ethanol Gas Standard Cylinder Report (SFN59282, Form 105-G).
 - B. ACA Test Record.

Chemical Test Operator Signature <i>Cory Mortensen</i>	Date <i>1-20-24</i>
Reviewed By (Crime Laboratory Use Only) <i>Aime Nahl</i>	Date 23Jan2024

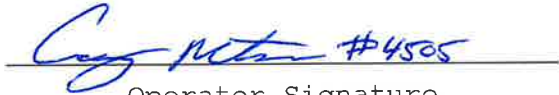
Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-006494
Location = CAVA 8164.14.00 09/16
01/20/2024 08:57

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	08:58
02 Std. Gas	0.083	08:59
03 Room Air	0.000	08:59
04 Std. Gas	0.083	09:00
05 Room Air	0.000	09:00
06 Std. Gas	0.083	09:00
07 Room Air	0.000	09:01

Lot No = 14323080A4
Cyl No = 6
Exp Date = 06/05/2025
County = 34 Oper No. = 020832


Operator Signature
CORY MORTENSEN

Remarks:
ACA Test

Form 106-I8000