



INTOXILYZER® 8000 INSTALLATION AND REPAIR CHECKOUT
 NORTH DAKOTA OFFICE OF ATTORNEY GENERAL
 CRIME LABORATORY DIVISION-TOXICOLOGY SECTION/BREATH ALCOHOL PROGRAM
 SFN 59281 (06/2018)

Serial Number 80-064191	Instrument Location ND SU PD-SBCP
Reason for Install/Repair <input checked="" type="checkbox"/> Install After Receiving From Crime Laboratory <input type="checkbox"/> Install After Location Change <input type="checkbox"/> Other (Specify) _____	

Check When Done:

- 1. Surge Protector Installed/Property Grounded.
- 2. Telephone Line Connected to Intoxilyzer® 8000.
- 3. Breath Tube Heated.
- 4. Enter Preliminary Data (i.e. Date, Time, DST (Y), and Location; Level 2, Function E).
- 5. Scan/Enter Gas Cylinder Information (Level 1, Function S).
- 6. Run Tests:
 - A. Print Test (Level 1, Function P).
 - B. ACA Test (Level 1, Function C).
 - C. Radio Frequency Interference (RFI) Test (CMS Mode or Level 1, Function B or C; Key Radio During Test).
- 7. Repair and/or Maintenance Performed (if any): _____
- 8. Complete the Top Portion of the Intoxilyzer® Record (SFN50496, Form 120-G) and Place it by the Intoxilyzer® for Use.
- 9. File Previous Intoxilyzer® Record (SFN504096, Form 120-G) at the Intoxilyzer® Location at the Agency.
- 10. Send the Following to the Crime Laboratory: Completed Intoxilyzer® 8000 Installation and Repair Checkout (SFN59281, Form 104-G), Print Test, ACA Test, and RFI Test.

Field Inspector Signature <i>Adam M. [Signature]</i>	Date 8/30/24
---	------------------------

Crime Laboratory Use Only

This installation has been reviewed and the instrument is approved to be used for the analysis of breath to determine alcohol concentration from the date the Field Inspector performed the installation. This record on file at the Office of Attorney General, Crime Laboratory Division, in the County of Burleigh, North Dakota, is certified to be a true and correct copy of the documents received.

Reviewed/Certified By <i>Anna Narehood</i>	Anna Narehood	Certified Date 03Sep2024
---	---------------	------------------------------------

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004191
Location = SBCP 8164.14.00 09/16
08/29/2024 18:47

***** Printer Test *****

abcdefghijklmnopqrstuvwxy1234567890- = |
ABCDEFGHIJKLMNopqrstuvwxyz!@#%&^* () _ + ?

abcdefghijklmnopqrstuvwxy1234567890- = |
ABCDEFGHIJKLMNopqrstuvwxyz!@#%&^* () _ + ?

Current Instrument Setup

Data Entry Mode: Enabled
Start Test Sequence: DABACABA
Display Prelim Rslt? Yes
Display Third Digit? Yes
Inhib Printer(Y/N)? No
Display Volume? No
Disable On Memfull? Yes
of Print Copies? 1
Select Std (D/W/I)? Dry

Standard Value? 0.075
Standard Lot #? 28423080A3
Standard Cyl #? 2
Standard Expiration? 11/05/2025
Oper No? 129972

Flow Cal. Date: 01/28/2010
Slope 738
Intercept -777574

IR Calibration Date: 06/25/2024
 3um 9um

0th Coef(*100): -15752 -25149
1st Coef(*100): 278540 137807
2nd Coef(*100): 2188 1022
H2O adj(mg/l*10k): 202 388

***** Printer Test End *****



Operator Signature
ADAM MALAFA

Remarks:

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501


CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004191
Location = SBCP 8164.14.00 09/16
08/29/2024 18:47

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	18:48
02 Std. Gas	0.079	18:48
03 Room Air	0.000	18:49
04 Std. Gas	0.079	18:49
05 Room Air	0.000	18:50
06 Std. Gas	0.078	18:50
07 Room Air	0.000	18:51

Lot No = 28423080A3
Cyl No = 2
Exp Date = 11/05/2025

County = 09 Oper No. = 129972



Operator Signature
ADAM MALAFA

Remarks: *Install Test*

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501


CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004191
Location = SBCP 8164.14.00 09/16
08/29/2024 18:51

DRY CAL CHECK

Test	AC	Time
01 Room Air	RFI*	18:52
02 Room Air	RFI*	18:52

*Invalid Test
Inhibited - RFI

Lot No = 28423080A3
Cyl No = 2
Exp Date = 11/05/2025
County = 09 Oper No. = 129972

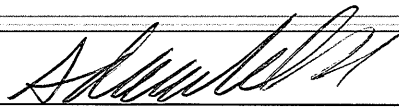

Operator Signature
ADAM MALAFA

Remarks: *Install RFI Test*

Form 106-I8000