



ETHANOL GAS STANDARD CYLINDER REPORT

NORTH DAKOTA OFFICE OF ATTORNEY GENERAL

CRIME LABORATORY DIVISION-TOXICOLOGY SECTION/BREATH ALCOHOL PROGRAM

SFN 59282 (08/2018)

Chemical Test Operator Name (Print) <i>Thomas, Patrick S.</i>		
Location <i>Fargo Police Department</i>	Intoxilyzer® Serial Number <i>80-004942</i>	
Gas Lot Number <i>01524080A4</i>	Gas Cylinder Number <i>32</i>	Gas Expiration Date <i>03-05-2026</i>

Check When Done:

- 1. Scan/Enter Gas Cylinder Information (Level 1, Function S).
- 2. Perform an ACA Test (Level 1, Function C).
Write in the ACA Test Results (Reported to 3 Digits, Ex. 0.081).
A. *0.082* AC
B. *0.080* AC
C. *0.081* AC
- 3. Complete the Top Portion of the Intoxilyzer® Record (SFN50496, Form 120-G) and Place it by the Intoxilyzer® for Use.
- 4. File Previous Intoxilyzer® Record (SFN50496, Form 120-G) at the Intoxilyzer® Location at the Agency.

- AEN 5. Send the Following to the Crime Laboratory:
- A. Ethanol Gas Standard Cylinder Report (SFN59282, Form 105-G).
 - B. ACA Test Record.

Chemical Test Operator Signature <i>P. Thomas</i>	Date <i>01/15/2025</i>
Reviewed By (Crime Laboratory Use Only) Anna Narehood <i>Anna Narehood</i>	Date 16Jan2025


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004942
Location = FAPD 8164.14.00 09/16
01/15/2025 19:25

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	19:25
02 Std. Gas	0.082	19:26
03 Room Air	0.000	19:26
04 Std. Gas	0.080	19:27
05 Room Air	0.000	19:27
06 Std. Gas	0.081	19:27
07 Room Air	0.000	19:28

Lot No = 01524080A4
Cyl No = 032
Exp Date = 03/05/2026
County = 09 Oper No. = 130056

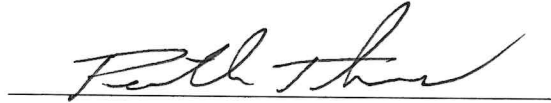

Operator Signature
PATRICK THOMAS

Remarks: *ACA for new Gas Cylinder*

Form 106-I8000