



ETHANOL GAS STANDARD CYLINDER REPORT

NORTH DAKOTA OFFICE OF ATTORNEY GENERAL

CRIME LABORATORY DIVISION-TOXICOLOGY SECTION/BREATH ALCOHOL PROGRAM

SFN 59282 (08/2018)

Chemical Test Operator Name (Print) <i>Tom Sharp</i>		
Location <i>MANDAN</i>	Intoxilyzer® Serial Number <i>80-006506</i>	
Gas Lot Number <i>14323080A1</i>	Gas Cylinder Number <i>73</i>	Gas Expiration Date <i>6/5/25</i>

Check When Done:

- ☒ 1. Scan/Enter Gas Cylinder Information (Level 1, Function S).
- ☒ 2. Perform an ACA Test (Level 1, Function C).
Write in the ACA Test Results (Reported to 3 Digits, Ex. 0.081).
A. 0.*081* AC
B. 0.*090* AC
C. 0.*081* AC
- ☒ 3. Complete the Top Portion of the Intoxilyzer® Record (SFN50496, Form 120-G) and Place it by the Intoxilyzer® for Use.
- ☒ 4. File Previous Intoxilyzer® Record (SFN50496, Form 120-G) at the Intoxilyzer® Location at the Agency.
- ☒ 5. Send the Following to the Crime Laboratory:
A. Ethanol Gas Standard Cylinder Report (SFN59282, Form 105-G).
B. ACA Test Record.

AEN

Chemical Test Operator Signature <i>[Signature]</i>	Date <i>4/3/25</i>
Reviewed By (Crime Laboratory Use Only) Anna Narehood <i>[Signature]</i>	Date 15April2025


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-006506
Location = MAND 8164.14.00 09/16
04/03/2025 14:26

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	14:27
02 Std. Gas	0.081	14:27
03 Room Air	0.000	14:28
04 Std. Gas	0.080	14:28
05 Room Air	0.000	14:29
06 Std. Gas	0.081	14:29
07 Room Air	0.000	14:30

Lot No = 14323080A1
Cyl No = 73
Exp Date = 06/05/2025
County = 30 Oper No. = 020611

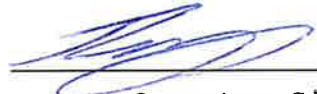

Operator Signature
TOM SHARP

Remarks:

Form 106-I8000