



# ETHANOL GAS STANDARD CYLINDER REPORT

NORTH DAKOTA OFFICE OF ATTORNEY GENERAL

CRIME LABORATORY DIVISION-TOXICOLOGY SECTION/BREATH ALCOHOL PROGRAM

SFN 59282 (08/2018)

Chemical Test Operator Name (Print) <u>Jordan Malchose</u>		
Location <u>Mott</u>	Intoxilyzer® Serial Number <u>80-007091</u>	
Gas Lot Number <u>08425080A1</u>	Gas Cylinder Number <u>4</u>	Gas Expiration Date <u>05/05/2027</u>

Check When Done:

☒ 1. Scan/Enter Gas Cylinder Information (Level 1, Function S).

☒ 2. Perform an ACA Test (Level 1, Function C).

Write in the ACA Test Results (Reported to 3 Digits, Ex. 0.081).

A. 0.079 AC

B. 0.079 AC

C. 0.080 AC

☒ 3. Complete the Top Portion of the Intoxilyzer® Record (SFN50496, Form 120-G) and Place it by the Intoxilyzer® for Use.

☒ 4. File Previous Intoxilyzer® Record (SFN50496, Form 120-G) at the Intoxilyzer® Location at the Agency.

☒ 5. Send the Following to the Crime Laboratory:

A. Ethanol Gas Standard Cylinder Report (SFN59282, Form 105-G).

B. ACA Test Record.

Chemical Test Operator Signature <u>[Signature]</u>	Date <u>05/17/2025</u>
Reviewed By (Crime Laboratory Use Only) <u>Anna Narehood</u> <u>[Signature]</u>	Date <u>22May2025</u>

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-007091  
Location = MOTT      8164.16.00 09/18  
05/17/2025      17:26

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	17:27
02 Std. Gas	0.079	17:27
03 Room Air	0.000	17:27
04 Std. Gas	0.079	17:28
05 Room Air	0.000	17:28
06 Std. Gas	0.080	17:29
07 Room Air	0.000	17:29

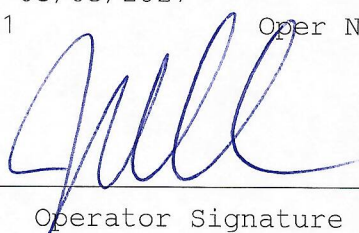
Lot No = 08425080A1

Cyl No = 04

Exp Date = 05/05/2027

County = 21

Oper No. = 132971

  
Operator Signature  
JORDAN MALCHOSE

Remarks:

Form 106-I8000