



ETHANOL GAS STANDARD CYLINDER REPORT

NORTH DAKOTA OFFICE OF ATTORNEY GENERAL

CRIME LABORATORY DIVISION-TOXICOLOGY SECTION/BREATH ALCOHOL PROGRAM

SFN 59282 (08/2018)

Chemical Test Operator Name (Print) <i>Taylor Wray</i>		
Location <i>BMD</i>	Intoxilyzer® Serial Number <i>80-005956</i>	
Gas Lot Number <i>19825080A2</i>	Gas Cylinder Number <i>13</i>	Gas Expiration Date <i>09/05/2027</i>

Check When Done:

- ☒ 1. Scan/Enter Gas Cylinder Information (Level 1, Function S).
- ☒ 2. Perform an ACA Test (Level 1, Function C).
Write in the ACA Test Results (Reported to 3 Digits, Ex. 0.081).
 - A. *0.016* AC
 - B. *0.016* AC
 - C. *0.016* AC
- ☒ 3. Complete the Top Portion of the Intoxilyzer® Record (SFN50496, Form 120-G) and Place it by the Intoxilyzer® for Use.
- ☒ 4. File Previous Intoxilyzer® Record (SFN50496, Form 120-G) at the Intoxilyzer® Location at the Agency.
- ☒ 5. Send the Following to the Crime Laboratory:
 - A. Ethanol Gas Standard Cylinder Report (SFN59282, Form 105-G).
 - B. ACA Test Record.

Chemical Test Operator Signature <i>Taylor Wray</i>	Date <i>09/09/2025</i>
Reviewed By (Crime Laboratory Use Only) <i>Anna Narehood</i>	Date <i>09/17/2025</i>

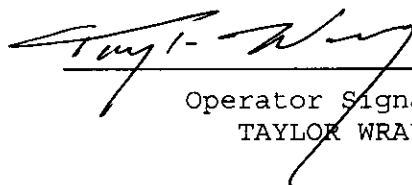
Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-005956
Location = BMDC 8164.14.00 09/16
09/09/2025 10:38

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	10:39
02 Std. Gas	0.076	10:39
03 Room Air	0.000	10:39
04 Std. Gas	0.076	10:40
05 Room Air	0.000	10:40
06 Std. Gas	0.076	10:41
07 Room Air	0.000	10:41

Lot No = 19825080A2
Cyl No = 13
Exp Date = 09/05/2027
County = 08 Oper No. = 131856


Operator Signature
TAYLOR WRAY

Remarks:

Form 106-I8000