

**INTOXILYZER® 8000 INSTALLATION AND REPAIR CHECKOUT**

NORTH DAKOTA OFFICE OF ATTORNEY GENERAL

CRIME LABORATORY DIVISION-TOXICOLOGY SECTION/BREATH ALCOHOL PROGRAM

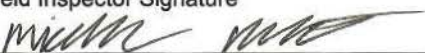
SFN 59281 (06/2018)

Serial Number 80-006490	Instrument Location ASHL
Reason for Install/Repair ☒ Install After Receiving From Crime Laboratory ☐ Install After Location Change ☐ Other (Specify) _____	

Check When Done:

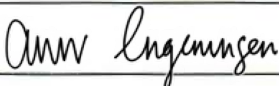
AEI
08Sept26

- ☒ 1. Surge Protector Installed/Property Grounded.
- ☒ 2. Telephone Line Connected to Intoxilyzer® 8000.
- ☒ 3. Breath Tube Heated.
- ☒ 4. Enter Preliminary Data (i.e. Date, Time, DST (Y), and Location; Level 2, Function E).
- ☒ 5. Scan/Enter Gas Cylinder Information (Level 1, Function S).
- ☒ 6. Run Tests:
 - ☒ A. Print Test (Level 1, Function P).
 - ☒ B. ACA Test (Level 1, Function C).
 - ☒ C. Radio Frequency Interference (RFI) Test (CMS Mode or Level 1, Function B or C; Key Radio During Test).
- ☐ 7. Repair and/or Maintenance Performed (if any): _____
- ☒ 8. Complete the Top Portion of the Intoxilyzer® Record (SFN50496, Form 120-G) and Place it by the Intoxilyzer® for Use.
- ☒ 9. File Previous Intoxilyzer® Record (SFN504096, Form 120-G) at the Intoxilyzer® Location at the Agency.
- ☒ 10. Send the Following to the Crime Laboratory: Completed Intoxilyzer® 8000 Installation and Repair Checkout (SFN59281, Form 104-G), Print Test, ACA Test, and RFI Test.

Field Inspector Signature 	Date 9-8-26
---	-----------------------

Crime Laboratory Use Only

This installation has been reviewed and the instrument is approved to be used for the analysis of breath to determine alcohol concentration from the date the Field Inspector performed the installation. This record on file at the Office of Attorney General, Crime Laboratory Division, in the County of Burleigh, North Dakota, is certified to be a true and correct copy of the documents received.

Reviewed/Certified By Anna Ingemansen 	Certified Date 08Sept2026
--	-------------------------------------

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-006490
Location = ASHL 8164.14.00 09/16
09/08/2026 10:52

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	10:54
02 Std. Gas	0.078	10:54
03 Room Air	0.000	10:54
04 Std. Gas	0.078	10:55
05 Room Air	0.000	10:55
06 Std. Gas	0.077	10:56
07 Room Air	0.000	10:56

Lot No = 264437
Cyl No = 111
Exp Date = 01/29/2028
County = 26 Oper No. = 131100



Operator Signature
MICHAEL MAROLT

Remarks:

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-006490
Location = ASHL 8164.14.00 09/16
09/08/2026 10:58

DRY CAL CHECK

Test	AC	Time
01 Room Air	RFI*	10:58
02 Room Air	0.000	10:59

*Invalid Test
Inhibited - RFI

Lot No = 264437
Cyl No = 111
Exp Date = 01/29/2028
County = 26 Oper No. = 131100



Operator Signature
MICHAEL MAROLT

Remarks:

Form 106-I8000