



ETHANOL GAS STANDARD CYLINDER REPORT

NORTH DAKOTA OFFICE OF ATTORNEY GENERAL

CRIME LABORATORY DIVISION-TOXICOLOGY SECTION/BREATH ALCOHOL PROGRAM

SFN 59282 (08/2018)

Chemical Test Operator Name (Print) <i>MICHAEL MAROLT</i>		
Location <i>NAPD</i>		Intoxilyzer® Serial Number <i>80-004200</i>
Gas Lot Number <i>264944</i>	Gas Cylinder Number <i>298</i>	Gas Expiration Date <i>5-12-2028</i>

Check When Done:

- AEI 21Sep2026
- ☒ 1. Scan/Enter Gas Cylinder Information (Level 1, Function S).
 - ☒ 2. Perform an ACA Test (Level 1, Function C).
Write in the ACA Test Results (Reported to 3 Digits, Ex. 0.081).

A. *0.081* AC
B. *0.080* AC
C. *0.081* AC

- ☒ 3. Complete the Top Portion of the Intoxilyzer® Record (SFN50496, Form 120-G) and Place it by the Intoxilyzer® for Use.
- ☒ 4. File Previous Intoxilyzer® Record (SFN50496, Form 120-G) at the Intoxilyzer® Location at the Agency.
- ☒ 5. Send the Following to the Crime Laboratory:
 - A. Ethanol Gas Standard Cylinder Report (SFN59282, Form 105-G).
 - B. ACA Test Record.

Chemical Test Operator Signature <i>[Signature]</i>	Date <i>09-21-2026</i>
Reviewed By (Crime Laboratory Use Only) Anna Ingemansen <i>Anna Ingemansen</i>	Date 21Sep2026

Checklist #2 was checked by Anna Ingemansen on 21Sep2026. The ACA test record was attached. This does not affect the Installation of the instrument. 21Sep2026, AEI


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004200
Location = NAPO 8164.14.00 09/16
09/21/2026 10:35

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	10:36
02 Std. Gas	0.081	10:36
03 Room Air	0.000	10:37
04 Std. Gas	0.080	10:37
05 Room Air	0.000	10:38
06 Std. Gas	0.081	10:38
07 Room Air	0.000	10:39

Lot No = 264944
Cyl No = 298
Exp Date = 05/12/2028
County = 24 Oper No. = 131100

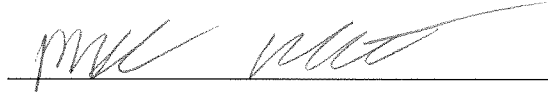

Operator Signature
MICHAEL MAROLT

Remarks:

Form 106-I8000